



EASY-Spray PepMap RSLC C18 Column

Increased sample throughput with a capillary EASY-Spray column

Benefits

- Higher sensitivity compared with Analytical scale LC
- Higher throughput than in standard nanoLC experiments
- Robust column design for high throughput proteomics
- Designed for TFA-free LC-MS analysis, minimizing ion suppression effects
- Easy to use format for Capillary LC

Keywords

Capillary LC, high throughput, proteomics

Introduction

The Thermo Scientific™ EASY-Spray™ 2 μm 15 cm \times 150 μm capillary column raises the bar in proteomics by greatly increasing the number of samples that can be analyzed by proteomics laboratories when compared to classical nanoLC methods. Capillary LC offers advantages in terms of ease of use when compared to nanoLC and also higher sensitivity compared to Analytical LC. The column is fully compatible with Thermo Scientific™ Capillary LC and Thermo Scientific™ Orbitrap™ mass spectrometer systems to offer a seamless integration of speed, sensitivity and robustness required for high throughput proteomics workflows.

The EASY-Spray design is a simple and easy to use “plug and spray” approach designed to improve throughput and is a novel alternative to classic low flow LC columns and sources. Integrating this technology with a 150 μm i.d. column further enhances productivity by the ability to increase throughput capacity for current nanoLC users, or by increasing sensitivity for analytical scale LC users.

Column technology

The EASY-Spray 15 cm x 150 μm column is packed with Thermo Scientific™ Acclaim™ PepMap™ C18 2 μm 100Å particles and has a 150 μm inner diameter. It is used for the separation of digested peptides that are used for the identification of proteins in proteomics applications. The Acclaim PepMap column is designed for TFA-free analysis, which minimizes ion suppression effects in mass spectrometry experiments.

Applications

Improved sensitivity

There is a clear experimental relationship between flow rate and sensitivity that can be described by a power law function (Figure 1). Thus, major improvements in sensitivity can be achieved at nano and capillary flow rates compared to analytical flow rates.

The EASY-Spray capillary column with 150 μm ID can be operated at a wide flow range from 1 to 3 $\mu\text{L}/\text{min}$, corresponding to linear velocities from 2 to 4 mm/s. Results have shown that capLC-MS can be 40 times more sensitive than analytical flow LC-MS at 450 $\mu\text{L}/\text{min}$ and between 2 to 4 times less sensitive than nanoLC-MS at 300 nL/min (Figure 1), when loading the same amount of sample onto the column*.

* Sensitivity gain or loss are dependent on compound and application

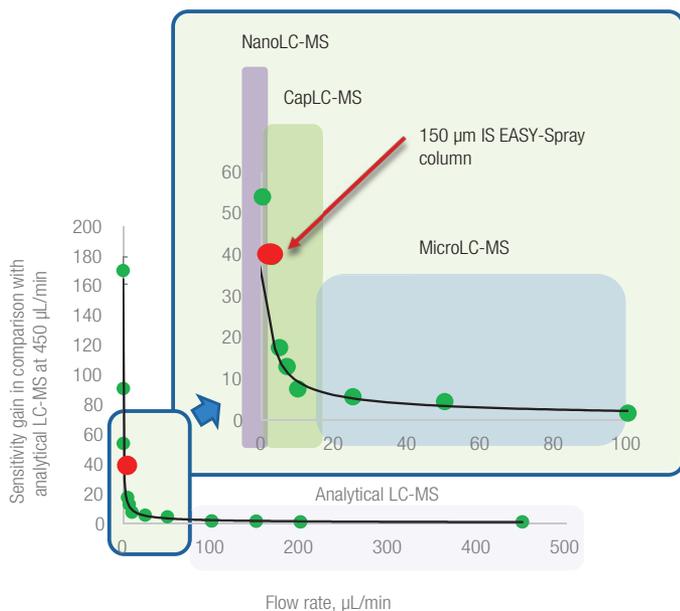


Figure 1. The increase in sensitivity experienced with low-flow LC-MS compared to analytical flow LC-MS

The EASY-Spray 15 cm x 150 μm column has an optimal flow rate of 3 $\mu\text{L}/\text{min}$ when used for high-throughput, targeted analysis. When used in profiling of proteomics samples with our capLC-MS platform, an optimal flow rate of 1.2 $\mu\text{L}/\text{min}$ was been shown to be ideal for the application used. Typical chromatograms for capLC-MS separations at 3 and 1.2 $\mu\text{L}/\text{min}$ (Figure 2 and 3) show that the separation performance is not compromised at low flow rates and is comparable to analytical flow separations. The peak width at half maximum (PWHM) is about 3 s for fast capLC-MS separations at 3 $\mu\text{L}/\text{min}$ and less than 15 s for longer gradients at 1.2 $\mu\text{L}/\text{min}$, which results in peak capacities of approx. 300 for both gradients.

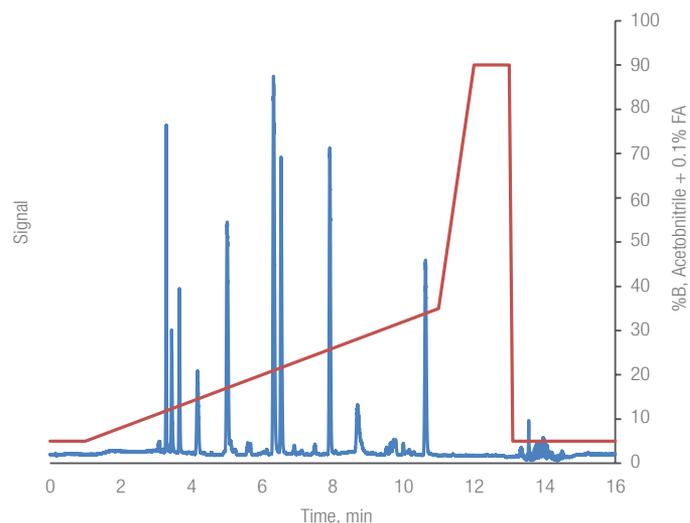


Figure 2. Typical chromatogram for 1 pmol of CytC protein digest on column. Flow rate: 3 $\mu\text{L}/\text{min}$ backpressure: 450 bar

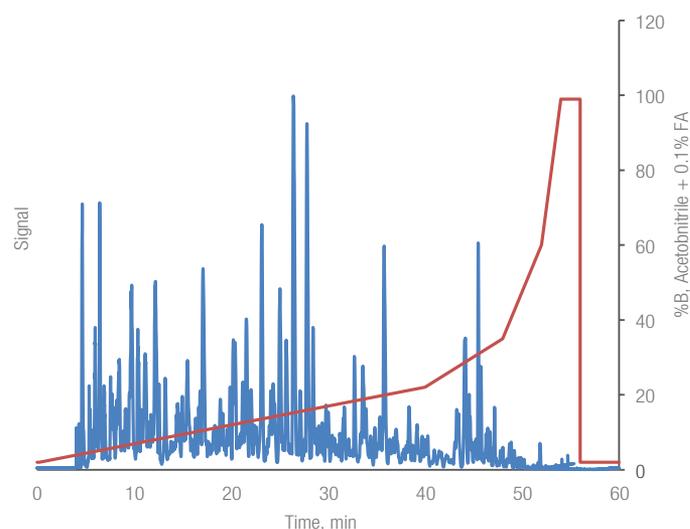


Figure 3. Typical chromatogram for 1 μg of HeLa cell lysate digest on column. Flow rate: 1.2 $\mu\text{L}/\text{min}$ backpressure: 200 bar

Robustness

The EASY-Spray 15 cm × 150 μm column was monitored over an extended injection sequence. A Thermo Scientific™ capillary-flow UltiMate™ 3000 RSLCnano system capLC system was configured for direct injections and evaluated over 8 days of operation (Figure 4a).

Over this time period, an excellent retention time stability (RSD < 1%) was observed for 350 injections.

The peak area stability was evaluated on MS1 level for 150 consecutive injections conducted between day 5 and 8 (Figure 4b). The peak areas achieved RSD values of less than 10% and were obtained even without internal standard correction during long term testing.

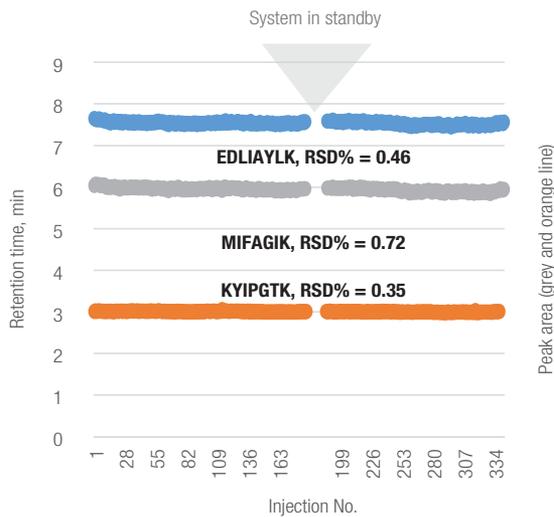


Figure 4a. Retention time stability for cap LC-MS analysis of CytC protein digest.

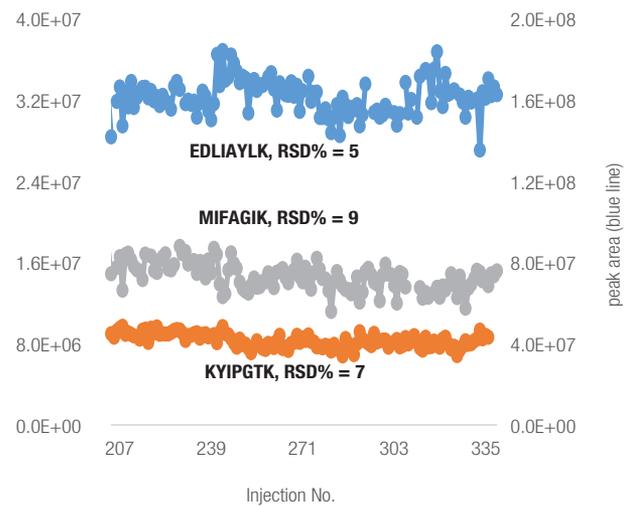


Figure 4b. Peak area stability for cap LC-MS analysis of CytC protein digest.

Operational specifications

Parameter	Recommendation
Typical applications	Peptide analyses
Base material	Spherical silica, 100 Å pores
Particle size	2 µm
Chemistry	C18, endcapped
Column dimensions	150 µm i.d. × 15 cm
Recommended flow	1–3 µL/min
Recommended temperature	25–60 °C
Maximum pressure	1000 bar
Recommended sample quantity ¹	5 µg
pH stability	2–8
Solvent compatibility	All common RP solvents

¹ The recommended sample quantity of the Acclaim PepMap100 150 µm i.d. columns is given in weight amount of a protein digest of bovine serum albumin.

Ordering Information

Description	Particle size (µm)	Part number
EASY-Spray PepMap RSLC C18 15 cm × 150 µm	2	ES906

Find out more at thermofisher.com/EASYSpray

For Research Use Only. Not for use in diagnostic procedures. ©2020 Thermo Fisher Scientific Inc.
All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries. This information is presented as an example of the capabilities of Thermo Fisher Scientific products. It is not intended to encourage use of these products in any manners that might infringe the intellectual property rights of others. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representatives for details. **PS21714-EN 0820**

ThermoFisher
SCIENTIFIC