



Biosafety cabinets

Protection and convenience

1500 Series A2
Biological Safety Cabinet

Sample protection and user safety

Containment, comfort and convenience



The Thermo Scientific™ 1500 Series Class II, Type A2 biological safety cabinets combine first-class protection, energy efficiency and user convenience for today's modern laboratory. The cabinet is designed to offer excellent personal and product protection from biological hazards and contamination. It can be connected to an external exhaust to allow work with volatile toxic chemicals that are sometimes used in the cabinet.

The touchscreen displays the health status of the cabinet, provides performance and background data and gives the possibility to program the UV light at the end of the work period.



Containment

- Thermo Scientific™ SmartFlow™ Plus technology automates airflow compensation
- Digital Airflow Verification (DAVe™ Plus) alerts user to out-of-tolerance airflows that could compromise product and operator protection
- H14 HEPA filters remove airborne contaminants
- Tactile window mechanism provides positive confirmation that the correct working height is selected for the sash window



Comfort

- Increased light intensity with LED lighting for a brighter work environment
- Low noise level enhances user focus
- 10° sloped front for exceptional comfort and ergonomics
- Spacious work area with comfortable armrests



Convenience


- 7" touch screen graphical user interface displays safety and performance data
- Thermo Scientific SmartClean front window design to improve access to the work area for cleaning
- Thermo Scientific™ SmartPort™ design provides contained passage of tubing and cables into the work area
- Programmable UV light enables timed sterilization run to be preset before work commences

7" Touchscreen for operation information

Product status, analysis and UV light programming



The unit is in working mode, with glass sash in its working position, and both the downflow and inflow air speeds at operational speeds. The white outline around the fan logo at the bottom means the blowers are all on.

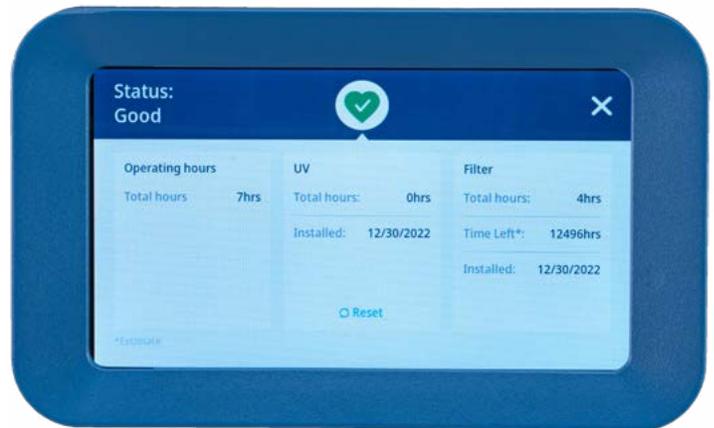

Tapping the 'Health Status Good' icon (the green heart) will tell you more the status of the unit – (1) the total accumulated operating hours of the unit, (2) the total UV work hours and the UV bulb install date, and (3) the filter running time, estimated residual lifetime of the HEPA filter, and the filter install date.



Clear visual and audible alarm are emitted from the device if the downflow, inflow, or other parameters are not at operational settings and operator is at risk of exposure to biohazardous materials.

Advanced features for user comfort



Energy efficiency

- DC motors offer low energy consumption
- Night Set-Back Mode reduces airflow to 30% during low-usage periods
- LED lighting
- Programmed UV auto-start and UV disinfection cycle duration

Specifications

Size / Width		4-foot	6-foot	
Dimensions	Exterior dimensions W x H x D	mm/in	1300 x 1568 x 800 / 51.2 x 61.8 x 31.5	1900 x 1568 x 800 / 74.8 x 61.8 x 31.5
	Interior dimensions W x H x D	mm/in	1200 x 780 x 630 / 47.2 x 30.7 x 24.8	1800 x 780 x 630 / 70.9 x 30.7 x 24.8
	Shipping dimensions ¹ W x H x D	mm/in	1410 x 1730 x 925 / 55.5 x 68.1 x 36.4	2010 x 1730 x 925 / 79.1 x 68.1 x 36.4
	Working height of front window	mm/in		250 / 10
	Maximum opening height of front window	mm/in		535 / 21
Weight	Net weight	kg/lbs	200 / 441	280 / 617
	Shipping weight ¹	kg/lbs	225 / 496	315 / 694
Electrical data	Voltage	V		230
	Frequency	Hz		50 / 60
Filter specification			H14 HEPA EN 1822, 99.995% @ MPPS, 99.999% at @ 0.3 µm	
Certification			NSF/ANSI 49, YY0569, UL, CE	
Ergonomics	Lighting power	W	20	40
	Average illumination	Lx	800 (SST) 900 (SmartCoat)	950 (SST) 1050 (SmartCoat)
	SmartPort access ²		Two 3" plugged cable ports, one on each side wall	
	Receptacles		Two standard single receptacles, located on rear wall, right and left sides	
	Service valves		Up to 6 total (3 on each side wall)	
10-inch front opening	Sound pressure level	dB (A)	<63	<65
	Energy consumption, operating set point ³	W	200	400
	Energy consumption, night set-back mode	W	70	120
	Heat output, operating mode (non-vented)	BTU/hr (W)	683 (200)	1366 (400)
	Heat output, night set-back mode (lights off, non-vented)	BTU/hr (W)	239 (70)	410 (120)
	Exhaust/inflow air volume	cfm (m ³ /h)	344 (585)	517 (878)
	Exhaust volume, thimble ducted	cfm (m ³ /h)	448 (761)	672 (1141)

¹ Export packaging dimensions and weight upon request

² Not available for catalog numbers: 1574 and 1579

³ Clean filters, fans operating at nominal set point, interior lighting activated

Ordering information

Bundle packages	Electrical configuration	Size	Certification	Cat. No. adjustable height stand	Cat. No. fixed height stand / castors
1500 Series Class II, Type A2 biological safety cabinet with one-piece stainless steel interior , UV light, armrests, and a stand.	China	4'	NSF & YY	1574-A	1574-AC
		6'	NSF & YY	1579-A	1579-AC
		4'	NSF	1574-A-GP	1574-AC-GP
		6'	NSF	1579-A-GP	1579-AC-GP
1500 Series Class II, Type A2 biological safety cabinet with SmartCoat interior , UV light, armrests, and a stand.		4'	NSF	1584-A	1584-AC
		6'	NSF	1589-A	1589-AC

You can also build your own 1500 Series A2 BSC, by starting with these cabinet-only part numbers:

Individual cabinets	Electrical configuration	Size	Certification	Cat. No. 10-inch front opening
1500 Series Class II, Type A2 biological safety cabinet with one-piece stainless steel interior	China	4'	NSF	1546
		6'	NSF	1548
1500 Series Class II, Type A2 biological safety cabinet with SmartCoat interior		4'	NSF	1556
		6'	NSF	1558

Packages are also available with inner receptacles and power cords that are commonly used in Asia:

Bundle packages	Electrical configuration	Size	Certification	Cat. No. adjustable height stand	Cat. No. fixed height stand / castors
1500 Series Class II, Type A2 biological safety cabinet with SmartCoat interior , UV light, armrests, and a stand.	Korea	4'	NSF	1584-F	1584-FC
		6'	NSF	1589-F	1589-FC
	Southeast Asia & Hong Kong	4'	NSF	1584-G	1584-GC
		6'	NSF	1589-G	1589-GC
	India	4'	NSF	1584-M	1584-MC
		6'	NSF	1589-M	1589-MC
1500 Series Class II, Type A2 biological safety cabinet with one-piece stainless steel interior , UV light, armrests, and a stand.	Korea	4'	NSF	1574-F	1574-FC
		6'	NSF	1579-F	1579-FC
	Southeast Asia & Hong Kong	4'	NSF	1574-G	1574-GC
		6'	NSF	1579-G	1579-GC
	India	4'	NSF	1574-M	1574-MC
		6'	NSF	1579-M	1579-MC

Factory-installed accessories

Accessory	Description	Cat. No.
UV light option	UV light option for 4 foot cabinet	141055-BSC
	UV light option for 6 foot cabinet (order 2x)	141055-BSC*

* Order two.

Field-installed accessories

Accessory	Description	Cat. No.
Stands	Manual adjustable stand for 4 foot cabinet, provides work surface height of 760 to 965 mm adjustable in 55 mm increments	1911302
	Manual adjustable stand for 6 foot cabinet, provides work surface height of 760 to 965 mm adjustable in 55 mm increments	1911303
	Electric adjustable stand for 4 foot cabinet (230V), work surface height adjustable from 750 mm to 950 mm	50109314
	Electric adjustable stand for 6 foot cabinet (230V), work surface height adjustable from 750 mm to 950 mm	50109325
	Rolling castor stand for 4 foot cabinet, provides fixed height of 32.5 inches	1921302
	Rolling castor stand for 6 foot cabinet, provides fixed height of 32.5 inches	1921303
Exhaust transitions	Thimble duct exhaust transition for 4 foot cabinet	1911316
	Thimble duct exhaust transition for 6 foot cabinet	1911317
	Alnor exhaust alarm for 1300 Series A2 biological safety cabinets, qualified install only	1910185
User convenience equipment	Interior tubing retainer clips (6 per kit)	1910083
	Hanging shelf for base stand	1910179
	Stainless steel armrest (set of two)	1911313
	Adjustable footrest	191127
	Replacement SmartPort grommets (4 per kit)	1910112
	Ergolign saddle stool, 23 inch to 30.5 inch adjustable seat height	ERG-CR-2330
	IV bag holder kit with 12 hooks (for 4 foot and 6 foot cabinet)	1911312
	Padded armrest 316 SS (set of 2 metal armrests plus gel inserts)	50156438
Replacement gel pads for padded armrest (2 pcs.)	50156437	
Alternate drain valve	Stainless steel ball valve (replaces standard non-corrosive PVC drain valve)	249098
Sidewall utilities	Combustible gas valve	1911308
	Non-combustible gas valve	1911309
	Vacuum valve	1911310
	Water valve	1911311
Rear cover kits*	Rear cover kits for 4 foot 1500 A2 cabinet	1915368
	Rear cover kits for 6 foot 1500 A2 cabinet	1915468

* Rear cover can also be factory-installed.



1910179 hanging shelf for base stand



ERG-CR-2330 Ergolign saddle stool



191127 optional adjustable footrest



Fixed height stand with castor



Stainless steel armrests



UV light



Learn more at thermofisher.com/1500bsc

For Laboratory Use. It is the customer's responsibility to ensure that the performance of the product is suitable for customers' specific uses or applications. © 2022-2023 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. EXT5605 0823