

Scale up

your biopharma
experiments

ThermoFisher
SCIENTIFIC



Orbitrap Ascend
BioPharma Tribrid
mass spectrometer



thermo scientific

Scale up biopharma insight

Set new standards in biotherapeutic characterization

Built to meet the demands of biopharmaceutical experiments, the Thermo Scientific™ Orbitrap™ Ascend BioPharma Tribrid™ mass spectrometer adds innovations to scale up characterization of biotherapeutic products and native proteins. Achieve greater coverage using a revolutionary hardware design that features two ion routing multipoles. Characterize intact biopharmaceuticals with improved Native MS workflows. Perform sensitive and deep peptide-based analyses and enhance elucidation of even the most difficult-to-characterize molecules, such as oligonucleotides, using a combination of optional orthogonal fragmentation techniques. With these capabilities and more, access unprecedented experimental throughput, versatility and usability to meet tomorrow's challenges while achieving more insight today.



“The Orbitrap Ascend BioPharma mass spectrometer combines two fundamental features that make it perfect for characterizing intact biopharmaceuticals. Multiple ion activation types are complemented by PTCR to simplify complex mass spectra and facilitate their interpretation. And the Native MS option increases the Orbitrap m/z range up to m/z 16,000. Together they enable researchers to perform sophisticated ion manipulations and analyze the most diverse molecules.”

Luca Fornelli, PhD
Professor of Biology
University of Oklahoma



Native proteomics

Characterize ADCs and other biotherapeutics



Peptide mapping

Perform PTM and disulfide mapping



Top-down characterization

Completely characterize protein complexes and proteoforms



Oligonucleotides


Fully characterize oligonucleotides



Orbitrap Ascend
BioPharma Tribrid mass spectrometer

Innovative technology addresses your biopharma needs

The Orbitrap Ascend BioPharma Tribrid mass spectrometer features innovations to scale up your biotherapeutic characterization with throughput, versatility and ease-of-use.

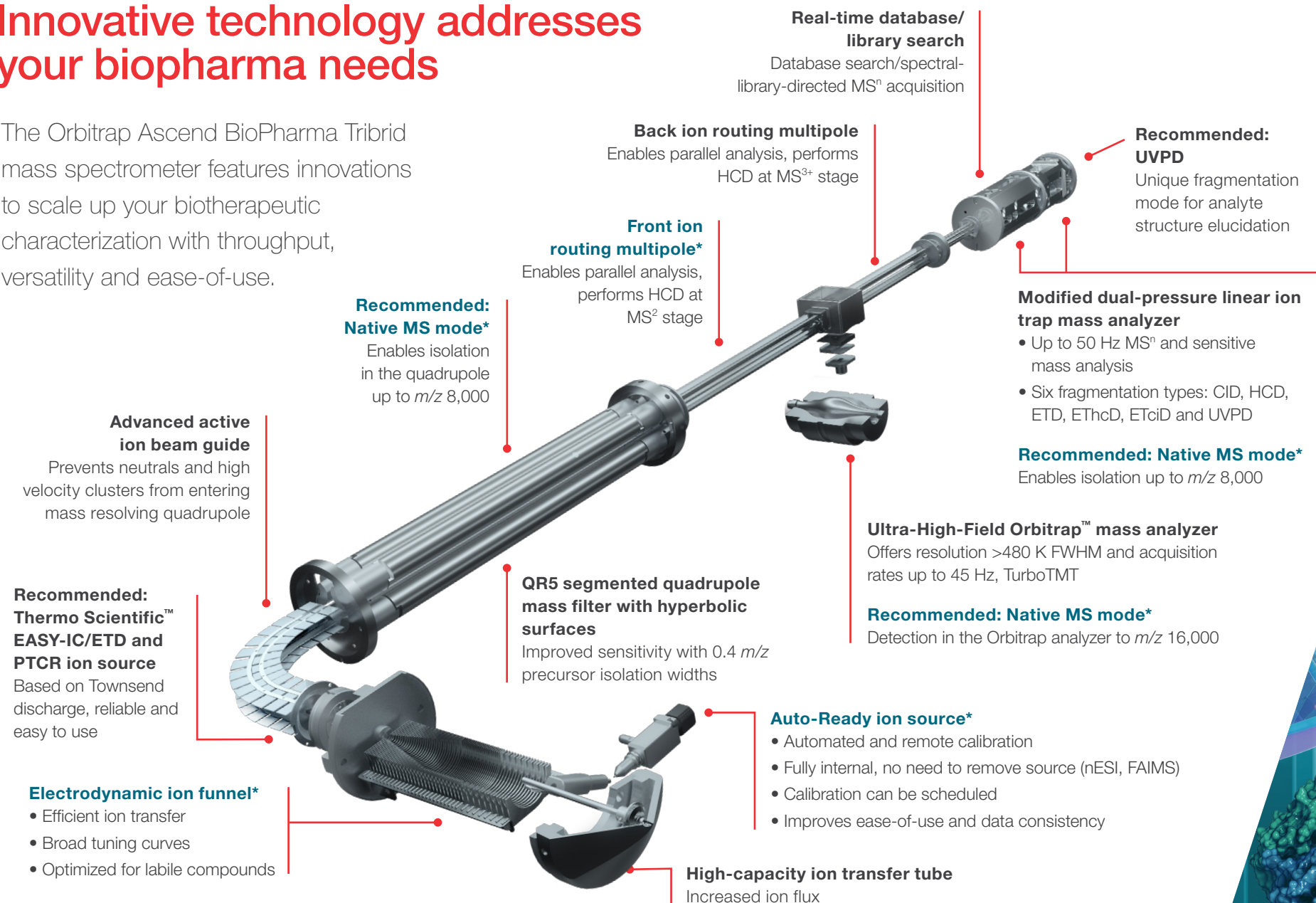

OPTIONS IC | ETD | PTCR | Native MS* | UVPD | FAIMS Pro Duo interface

*New on this platform

Peptides

Confidently map PTMs and characterize HCPs at analytical flow rates


Scale up to quickly characterize low abundance modifications

High mass and structural complexity make biotherapeutics susceptible to post-translational modifications (PTM), some of which can be detrimental to their safety, stability and efficacy. A versatile approach to PTM characterization, peptide mapping often requires the analysis of numerous samples, which challenges the analytical method to generate high-confidence data over a wide dynamic range in short liquid chromatography (LC) runtimes. Identifying deamidation is particularly difficult due to the small mass shift and the overlapping isotopic envelopes of unmodified and deamidated peptides.

With enhanced speed and sensitivity, the Orbitrap Ascend BioPharma Tribrid mass spectrometer can achieve 100% sequence coverage with identification of low abundance modifications—including deamidation—at analytical flow rates, increasing both sample throughput and certainty in results.

Orbitrap Ascend BioPharma Tribrid mass spectrometer analysis of Trastuzumab peptides at a flow rate of 300 μ L/min. A sequence coverage of 100% for the light and heavy chains of Trastuzumab at analytical flow rates was achieved. Low abundance modifications (<1% of the unmodified peptide's relative intensity) and the deamidation at 0.9% relative abundance were easily detected and characterized.


Proteins	Number of MS peaks	MS peak area	Sequence coverage
1:Heavy chain	2,420	46.3%	100.0%
2:Light chain	990	19.6%	100.0%

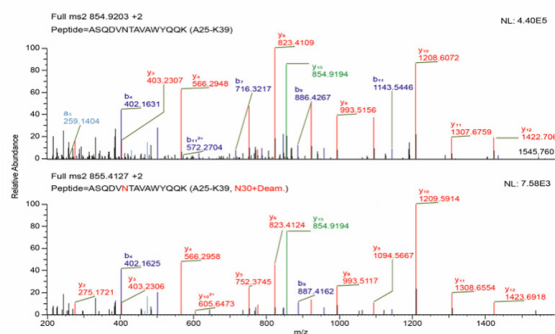

Scale up to confidently identify low-abundance HCPs

Host cell proteins (HCPs)—the biologic product impurities released during cell growth and processing—can detrimentally affect product safety and efficacy. LC with high-resolution accurate-mass (HRAM) MS analysis allows identification of multiple HCPs in a run and overcomes a typically huge intrasample dynamic range to detect very low abundance residual HCPs.

The sensitivity and dynamic range of the Orbitrap Ascend BioPharma Tribrid mass spectrometer provides rapid, high-confidence detection and identification of HCPs at less than 1 ppm concentration using analytical flow rates. Fragments generated via higher-energy collisional dissociation (HCD) enhance HCP identification.

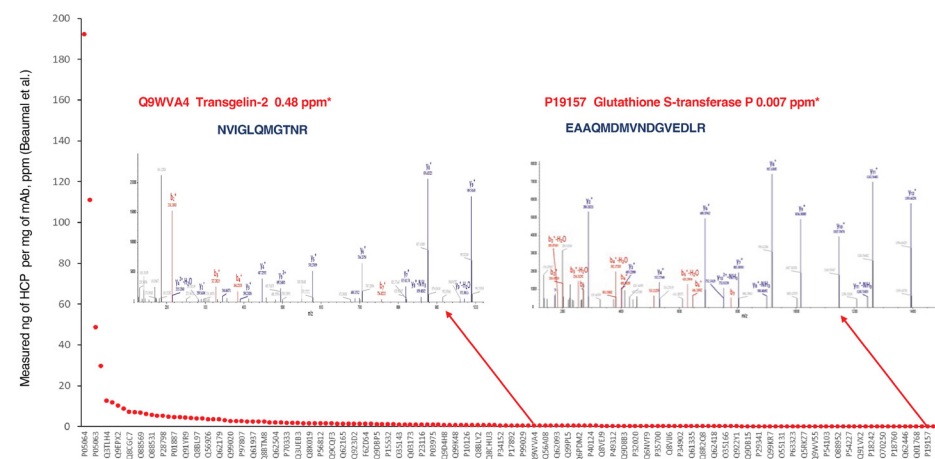

The Orbitrap Ascend BioPharma Tribrid mass spectrometer identified over 260 HCP peptides down to 7 ppb concentrations at 300 μ L/min flow rates.¹


1. Beaumont C, Beck A, Hernandez-Alba O et al. (2023) Advanced mass spectrometry workflows for accurate quantification of trace-level host cell proteins in drug products: Benefits of FAIMS separation and gas-phase fractionation DIA. *Proteomics and Systems Biology* 23 (16). doi.org/10.1002/psm.202300172.

ADC

Scale up native protein characterization

Protein therapeutics often generate undecipherable spectra due to their high levels of intrinsic heterogeneity. By using Orbitrap Ascend BioPharma Tribid mass spectrometer's recommended Native MS and proton transfer charge reduction (PTCR) options, interpretable spectra can be generated from heterogeneous protein samples. The Native MS option extends the quadrupole isolation range to m/z 2,000–8,000, which can be used in combination with narrow data-independent acquisition (DIA) windows to simplify the ion population in each scan. The Native MS option also extends Orbitrap analyzer detection to m/z 16,000. These capabilities enable unambiguous identification of previously indiscernible or low-abundance proteoforms. PTCR generates perfluoroperhydrophenanthrene (PFPP) ions for subsequent gas-phase, ion-ion reactions that produce lower charge state distributions to enhance protein characterization.

The flexibility of the Orbitrap Ascend BioPharma Tribid mass spectrometer enables scientists to carry out unique experiments that combine the recommended Native MS option with PTCR, as well as other dissociation techniques, making it a powerful instrument for comprehensive characterization of therapeutic proteins in their native state.

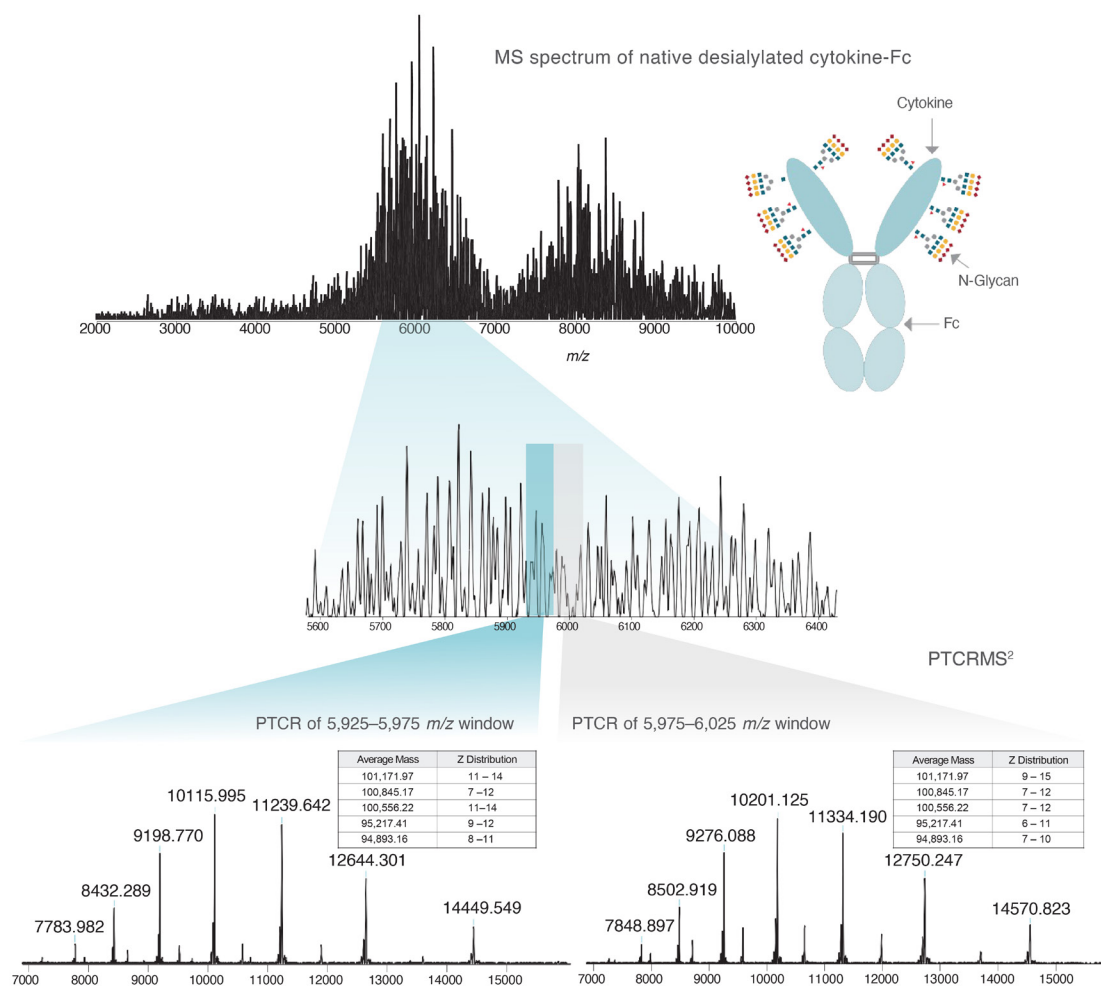

Orbitrap Ascend BioPharma Tribid mass spectrometer analysis of native desialylated cytokine-Fc analyzed with Native MS option. Two DIA windows were fragmented for comprehensive characterization.


Data courtesy of Wendy Sandoval.

Detect and identify unexpected low-abundance tRNA modifications

Oligos


Scale up oligonucleotide characterization

Therapeutic oligonucleotides produced by chemical synthesis can carry various types of impurities, some of which are present at very low levels. One method to achieve good separation and identification of modified intact tRNA is ion-pairing reversed-phase chromatography (IP-RPLC) coupled to HRAM MS analysis using the Orbitrap Ascend BioPharma Tribrid mass spectrometer. Intact detection of low-abundance oligonucleotide modifications can be paired with well-established digestion techniques for in-depth characterization using ion trap MS³-based collision-induced dissociation (CID) and the added versatility provided by recommended UVPD.

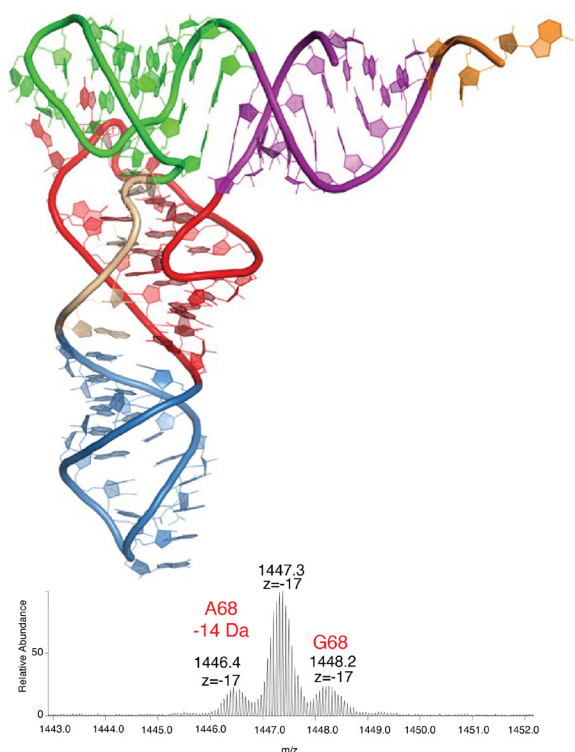


IP-RP coupled to Orbitrap Ascend BioPharma Tribrid mass spectrometer analysis enables separation of the modified tRNA species and identification of the known isodecoder and an under-modified variant.

Alternative fragmentation: MS³ by Resonance Activated CID in the Linear Ion Trap

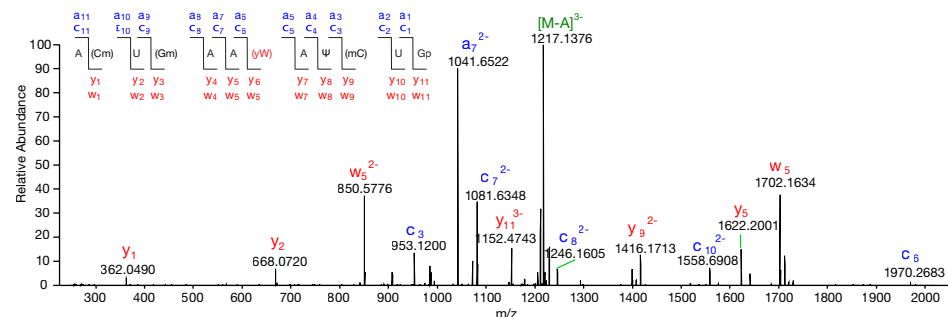

Alternative fragmentation: UVPD

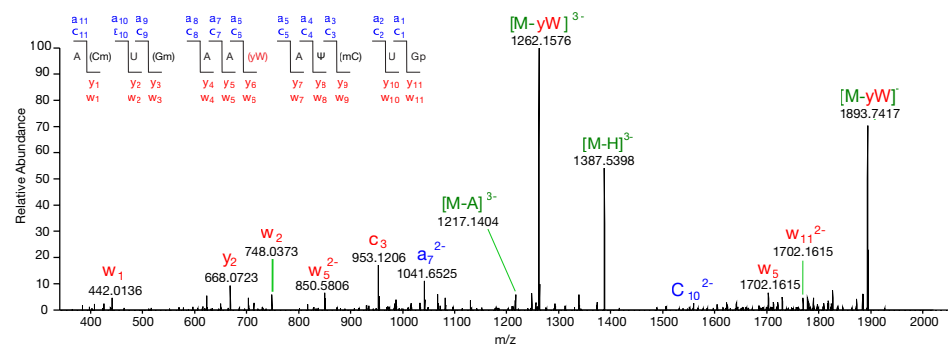

After digestion, the tRNA modifications within the oligonucleotide were located and well characterized. However, when subjected to CID, bulky side chains can be lost resulting in a large single peak limiting characterization of the oligonucleotide. Using alternative fragmentation, like recommended UVPD or MS³ targeting the base loss peak, results in an information-rich spectra.


Top-down

Achieve nearly 100% base-pair coverage for 100-mer sgRNA using low-q CID

Scale up characterization of sgRNA

Single guide RNA (sgRNA) plays a crucial role in the CRISPR-Cas9 gene editing system by guiding the Cas9 enzyme to specific genomic locations, enabling precise and targeted modifications of the DNA sequence. Characterization of sgRNA is essential to ensure precision from initial experiments to future therapies. Using the low frequency (q) CID capability of the Orbitrap Ascend BioPharma Tribrid mass spectrometer combined with intact and top-down methods, *de novo*-level fragmentation patterns of intact sgRNA can be obtained, providing nearly 100% sequence coverage for a 100-base-pair structure.

The low-q CID technique adds experimental flexibility by reducing precursor ion activation; the resulting product ion spectra are predominantly composed of sequence ions and base losses from the precursor. The reduction in product ion spectral complexity in favor of sequence-informative species simplifies and increases the confidence of m/z peak annotation and therefore confirmation or elucidation of the primary RNA sequence.

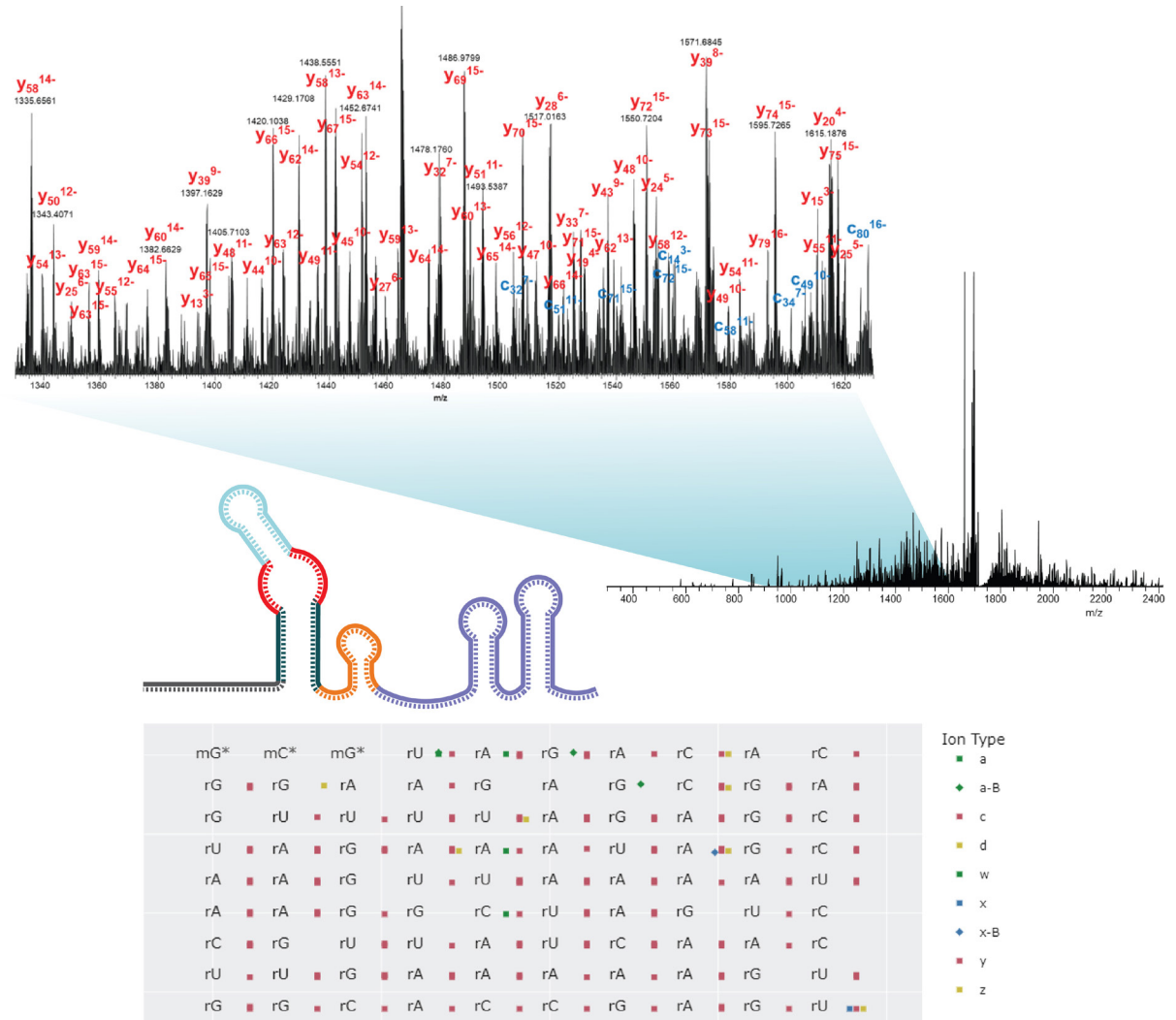

Orbitrap Ascend BioPharma Tribrid mass spectrometer analysis of a 100-mer sgRNA using low-q CID provided almost 100% sequence coverage.

Experience more high-quality results with less hassle using automated, remote and schedulable system checks and calibrations


Scale up convenience and ease-of-use

The Auto-Ready ion source is a fully integrated, standard, easy-to-use feature that increases laboratory productivity with automated, remote and schedulable system checks and internal calibrations. Because there is no need to remove the source (HESI, nESI or high-field asymmetric waveform ion mobility spectrometry [FAIMS]), there are no experimental setup interruptions required to perform internal calibrations. The user can automate the calibration to start at a scheduled time—for example, every week—when there are no experiments planned to run on the instrument. The calibration can run completely remotely, regardless of the nature of the last experiment. Because the calibration can be scheduled to occur regularly and automatically without interrupting vital work, users can expect to maintain mass spectrometer performance, improve data consistency and achieve more accurate and precise quantitation.




Auto-Ready ion source

Separate ion transfer tube

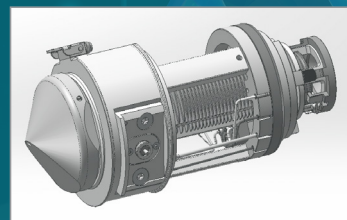



Dedicated emitter

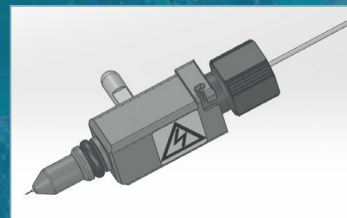

Robust delivery system



Automated weekly calibrations

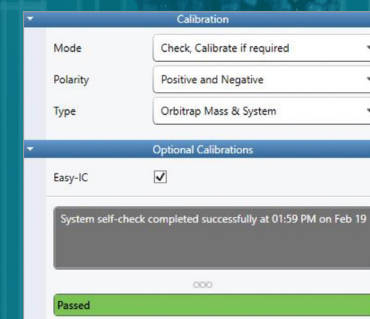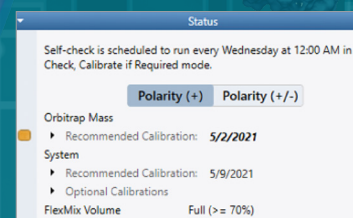

Resources and support



Services Central—All your service information at your fingertips

Spend less time searching for support and more time focusing on your important work. This online platform has what you need to easily manage your instruments and equipment.

Learn more at thermofisher.com/servicescentral



Technical and online support

Helping you keep your instruments running at peak performance is our goal. Whether you're looking for an instrument manual or spare parts, want to submit a repair request or check on the status of your warranty or service contract, we have every support option you're looking for.

Learn more at thermofisher.com/technicalresources



Protect your investments with expert lab services

Unity™ Lab Services provides a single source for integrated lab service, support and supply management. Our customized service offerings and world-class service experts have the flexibility and experience to address your laboratory's needs. We provide a complete portfolio of services and support solutions designed to help you improve productivity, reduce total cost of ownership and help ensure performance throughout your laboratory.

Learn more at unitylabservices.com

More resources

Thermo Scientific™ Ardia™ Platform
thermofisher.com/ArdiaProteomics

Thermo Scientific™ BioPharma Finder™ software
thermofisher.com/biopharmafinder

Thermo Scientific Compound Discoverer software
thermofisher.com/CompoundDiscoverer

mzCloud mass spectral library
thermofisher.com/mzCloud

Thermo Scientific Proteome Discoverer software
thermofisher.com/ProteomeDiscoverer

 Learn more at thermofisher.com/OrbitrapAscendBioPharma

General Laboratory Equipment – Not For Diagnostic Procedures. © 2024 Thermo Fisher Scientific Inc. All rights reserved. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. This information is presented as an example of the capabilities of Thermo Fisher Scientific products. It is not intended to encourage use of these products in any manner that might infringe the intellectual property rights of others. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details. **24-085-1011 | BR002863-EN 0524**

thermo scientific