


FreeZone® Tray Drying Accessories

Triad™, Stoppering, Bulk & Small Tray Dryers



726500—6-Port Tray Dryer Manifold. For use with FreeZone Stoppering and Bulk Tray Dryers. 30.0" w x 16.5" d x 7.0" h (76.2 x 41.9 x 17.8 cm). Six valves accommodate either 1/2" or 3/4" flask adapters. The 6-Port Tray Dryer Manifold and Isolation Valve 7761500 may not be field installed together. Shipping weight 19 lbs. (9 kg)

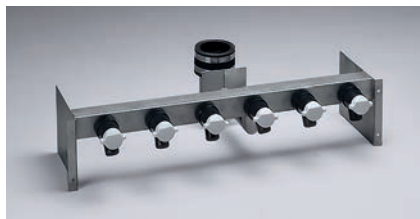


7204300—6-Port Tray Dryer Manifold and Isolation Valve. For use with FreeZone Stoppering and Bulk Tray Dryers. 30.1" w x 20.6" d x 13.5" h (76.5 x 52.3 x 34.3 cm). Six valves accommodate either 1/2" or 3/4" flask adapters. Shipping weight 45 lbs. (21 kg)


7761500—Isolation Valve. For use with FreeZone Stoppering and Bulk Tray Dryers. Separates tray dryer from freeze dryer to shorten time necessary to reestablish working vacuum levels and provide a means for checking end point. Consists of a valve, valve handle, clamps and couplings. The Isolation Valve and 6-Port Tray Dryer Manifold 726500 may not be field installed together. Shipping weight 5.0 lbs. (2.3 kg)

7756300—Shelf Spacers. For use with FreeZone Stoppering Tray Dryers. 9.0" w x 13.0" d x 2.5" h (22.9 x 33.0 x 6.4 cm). Type 304 stainless steel spacers placed on one or two empty shelves assist stoppering of partial loads of small serum bottles. One pair. Shipping weight 3.0 lbs. (1.4 kg)


Bulk Trays. Contain pre-frozen liquids for bulk drying. Type 304 stainless steel sheet with 0.9" (2.3 cm) lip. Shipping weight 3.0 lbs. (1.4 kg)

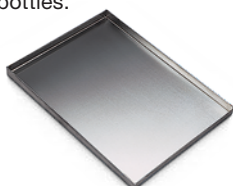

Catalog Number	Description	For Use With	Dimensions (w x d)
7377900	Small Tray	Stoppering, Bulk, Small	10.3" x 11.6" (26.2 x 29.5 cm)
7756200	Large Tray	Stoppering, Bulk	11.5" x 15.0" (29.2 x 38.1 cm)

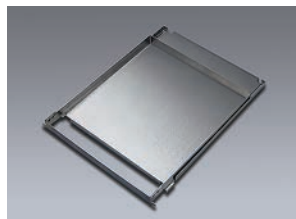


7756100—Tray with Slide-Out Bottom. For use with FreeZone Stoppering and Bulk Tray Dryers. 12.0" w x 14.0" d (30.5 x 35.6 cm). Stainless steel tray has separate bottom that slides out to allow glassware containers direct contact with shelf. Shipping weight 5.0 lbs. (2.3 kg)

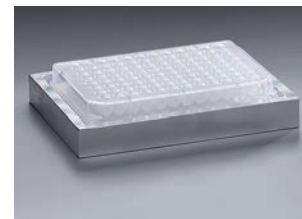

7369600—Microwell Plate Holder. Solid aluminum base holds a 127.8 x 85.5 mm microplate for excellent heat transfer and temperature stability for small samples susceptible to meltback. **Microwell plate not included.** 5.7" w x 8.0" d x 0.7" h (14.5 x 20.3 x 1.8 cm). One each. Shipping weight 2.0 lbs. (0.9 kg)

Rectangular Test Tube & Test Tube/Microcentrifuge Racks.


For use with any FreeZone Tray Dryer. Type 304 stainless steel rack, 12.0" w x 10.0" d x 2.5" h (30.5 x 25.4 x 6.4 cm), holds up to 70 each 12-13 mm test tubes. Test Tube/Microcentrifuge Rack has additional adapter plate for microcentrifuge tubes. **Tubes not included.**

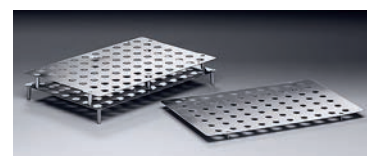

Catalog Number	Description	Shipping Weight
7370000	Test Tube Rack	2.8 lbs. (1.3 kg)
7370001	Test Tube/Microcentrifuge Tube Rack	3.5 lbs. (1.6 kg)



7439300—Microcentrifuge Tube Holder. 3.0" w x 2.25" d x 1.0" h (7.6 x 5.7 x 2.5 cm). Anodized aluminum block accommodates twelve 1.7 ml microcentrifuge tubes. Holder may be placed on any flat freeze drying surface. One each. **Microcentrifuge tubes are not included.** Shipping weight 0.6 lb. (0.3 kg)

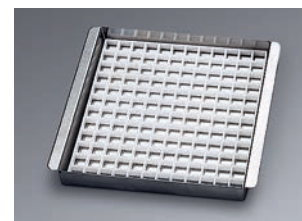

7516200—Support Grid. For use with any FreeZone Tray Dryer. 7.0" w x 7.0" d (17.8 x 17.8 cm). Stainless steel tray with removable plastic grid to hold ampules and other small specimen containers. Grid holds 144 containers of 12 mm dia. Shipping weight 2.5 lbs. (1.1 kg)

Authorised distributor

In Australia:

For customer service, call 1300-735-292

To fax an order, use 1800-067-639

To email an order, ordersau@thermofisher.com

In New Zealand:

For customer service, call 0800-933-966

To fax an order, use 0800-329-246

To email an order, ordersnz@thermofisher.com

ThermoFisher
S C I E N T I F I C