

# Cleanroom Consumables

Pharmaceutical | Biologics | Medical Device

Cleanroom Wipers | Sterile Products | Disinfectants | Mops | Swabs | Cleanroom Paper

## Table of Contents

---

Wipers.....	2-9
Vertex Wipers.....	2
Vertex Sterile Wipers .....	3
Sterile Wipers.....	4-5
Sterile Pre-Wetted Wipers.....	6-7
Dry Cleanroom Wipers .....	8-9
Pre-Wetted Cleanroom Wipers.....	10
Vertex Pre-Wetted Wipers.....	11
Sterile IPA.....	12
TexQ Disinfectant.....	13
Cleanroom Mops.....	14-15
Cleanroom and Cleaning Validation Swabs.....	16-17
TOC Cleaning Validation Kits .....	18
Cleanroom Labels.....	19
Cleanroom Paper .....	20-21
Cleanroom Tapes .....	22-23
Adhesive Mats .....	24
Product Index.....	25



### Product & Customer Support

Texwipe products are available through a worldwide network of distributors and are backed by our specialized technical support. We provide on-site facility audits, training posters, technical article reprints, product line brochures, technical datasheets, samples, test methods, and a comprehensive website.

## Critical Cleaning for Life Sciences

The Texwipe® team understands the unique requirements for cleaning and contamination control in the life sciences industries. These requirements extend beyond cleaning to cleaning validation and to the protocols to avert contamination.

Texwipe has a 40-year foundation of supplying the finest **cleanroom and sterile wipers** to life science companies. We are also the global leader in **sampling swabs** for cleaning validation and the only supplier of bound, **laboratory notebooks** designed for cleanrooms.

From **disinfectants** to **cleanroom tapes** to **Sterile IPA**, Texwipe offers everything Life Science companies need to ensure contamination control.

Texwipe products for pharmaceutical and biotechnology companies include:

- Wipers: Sterile, Dry, Pre-Wetted
- Sterile IPA
- TOC Cleaning Validation Kits
- Cleaning Validation Swabs
- Disinfectants
- Cleanroom Mops
- Adhesive Mats
- Cleanroom Paper
- Cleanroom Labels

Featured Product

### It's all about the numbers with disinfectants.

# TexQ™

Page 13

1-Step Cleaner & Disinfectant  
EPA Registered  
(10324-85)

2

EPA Registered  
(10324-85)

1

EPA Registered  
(10324-85)

1-Step Cleaner & Disinfectant  
EPA Registered  
(10324-85)

25%

more kill  
Claims

25%

more kill  
Claims

More detailed information about these or any other Texwipe products can be obtained by contacting your Texwipe representative at 1-336-996-7046 or by visiting us online at [www.texwipe.com](http://www.texwipe.com)



## What is Vertex®?

People are the greatest source of contamination in the cleanroom. We don't do this intentionally, it's simply a byproduct of the environment in which we live and our bodies sloughing off cells. By removing people from the manufacturing process, Texwipe has engineered a new level of clean which can't be touched... by us or the competition.

### No People.

From start to finish, this proprietary, patent-pending system isolates the wiper material from human contact. No people, no hands, no added risk of contamination.

### Consistency.

Vertex guarantees that we produce highly consistent wipers time after time. No variation in particulates or fibers, just the highest level of assurance that you will receive a highly consistent wiper, time after time.

### Options.

Vertex wipers were specifically developed with an ability to perform over a wide range of applications and settings. We offer a variety of fabrics in both 9" x 9" and 12" x 12" sizes, along with a selection Pre-Wetted Wipers.

### Vertex® Family of Dry Wipers

Product Number	Description	Packaging
TX22	High Durability, 12" x 12" nominal (30 cm x 30 cm), laser sealed-edge wipers	100 wipers/bag, 10 bags/case
TX29	High Durability, 9" x 9" nominal (23 cm x 23 cm), laser sealed-edge wipers	150 wipers/bag, 10 bags/case
TX42	High Sorption, 12" x 12" nominal (30 cm x 30 cm), laser sealed-edge wipers	100 wipers/bag, 10 bags/case
TX49	High Sorption, 9" x 9" nominal (23 cm x 23 cm), laser sealed-edge wipers	150 wipers/bag, 10 bags/case
TX59	Microdenier, 9" x 9" nominal (23 cm x 23 cm), sealed-edge wipers	100 wipers/bag, 10 bags/case

THE NEXT LEVEL OF CLEAN™

# Vertex® Sterile Wipers



STERILE

## A wiper so revolutionary, it's in a category by itself

Sterile Vertex Wipers are manufactured HandsFree™, eliminating the primary source of contamination risk.

### Lowest Risk of Microbial Contamination

By removing the human element, Vertex has the lowest risk for microbial contamination.

### Consistent Clean™

Minimum variability in particulates and fibers with our highly reproducible Roll2Bag™ process. Vertex offers the highest level of assurance that each wiper performs the same, bag after bag, lot after lot.

### High Sorption

100% polyester double-knit designed to provide a very soft feel plus high sorptive capacity.

Product Number	Description	Packaging
<b>TX3042</b>	Vertex Sterile High Sorption 12" x 12" sealed-edge wipers	100 wipers/bag (5 inner bags of 20 wipers), 5 bags/case
<b>TX3049</b>	Vertex Sterile High Sorption 9" x 9" sealed-edge wipers	100 wipers/bag (5 inner bags of 20 wipers), 5 bags/case

### Microdenier

Microdenier yarn has four times the number of filaments and twice the surface area in the yarn bundle as compared to standard polyester wipers. This increases wiping efficiency and creates the soft hand that is so uniquely characteristic of this microdenier wiper.

Product Number	Description	Packaging
<b>TX3052</b>	Vertex Sterile Microdenier 12" x 12" sealed-edge wipers	100 wipers/bag (5 inner bags of 20 wipers), 5 bags/case
<b>TX3059</b>	Vertex Sterile Microdenier 9" x 9" sealed-edge wipers	100 wipers/bag (5 inner bags of 20 wipers), 5 bags/case

### Pre-Wetted Wipers

USP-Purified IPA. Pre-wetted to just the right level, suitable for most applications.

Product Number	Description	Packaging
<b>TX3042P</b>	Sterile Vertex 12" x 12" sealed-edge pre-wetted with 70% USP Purified IPA / 30% DIW wipers	25 wipers/bag, 5 bags/case
<b>TX3044P</b>	Sterile Vertex® 12" x 12" sealed-edge wipers pre-wetted with 70% USP-Grade Denatured Ethyl Alcohol and USP Purified Water	25 wipers/bag, 5 bags/case
<b>TX3049P</b>	Sterile Vertex 9" x 9" sealed-edge pre-wetted with 70% USP Purified IPA / 30% DIW wipers	25 wipers/bag, 5 bags/case



# Sterile Dry Wipers

STERILE



## Aseptic and Sterile Areas Cleaned

Texwipe sterile wipers clean aseptic and sterile work areas efficiently while maintaining the integrity of the environments. They are ideal for cleaning aseptic fill areas, biotechnology production equipment, microbiological labs and isolators.

All Texwipe SterileWipe™ Wipers feature:

- **Gamma-irradiated to a 10<sup>-6</sup> sterility assurance level**
- **Specially designed packaging for aseptic environments**
- **Lot-specific Certificates of Irradiation and Compliance**
- **Clearly visible lot code and expiration date**
- **Protective strip on packaging to cover the wipers after opening**



## ISO Class 3-5

**SterileWipe™ LP 10**

Polyester sealed-border wiper with optimum balance of wiping efficiency, abrasion resistance and sorption.

Product Number	Description	Packaging
TX3212	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 5 inner bags of 20 wipers, 5 bags/case

**SterileWipe™ AS 10**

Premium two-ply, high-sorbency sealed-border wiper. Best choice for all standard applications.

Product Number	Description	Packaging
TX3215	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 5 inner bags of 20 wipers, 5 bags/case

**Sterile TexTra™ 10**

Two-ply capacity sorption in a single-ply, heavyweight polyester sealed-border wiper. All purpose wiper, excellent for spill control.

Product Number	Description	Packaging
TX3224	9" x 9" nominal (23 cm x 23 cm)	100 wipers/bag, 5 inner bags of 20 wipers, 5 bags/case
TX3225	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 5 inner bags of 20 wipers, 5 bags/case

## ISO Class 5-6

**Sterile TexTra™**

Two-ply capacity sorption in a single-ply, heavyweight, standard-edge polyester wiper. Economical and excellent for spill control.

Product Number	Description	Packaging
TX3220	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 5 inner bags of 20 wipers, 5 bags/case

**SterileWipe™ LP**

Standard-edge polyester all-purpose wiper. Combines low particle and extractable levels with superior strength.

Product Number	Description	Packaging
TX3211	9" x 9" nominal (23 cm x 23 cm)	100 wipers/bag, 10 inner bags of 20 wipers, 5 bags/case

## ISO Class 6-7

**SterileWipe™ HS II**

Absorbent, nonwoven cellulose/polyester wiper provides the optimum balance of sorption, durability and cleanliness. Economical all-purpose wiper.

Product Number	Description	Packaging
TX3210	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 5 inner bags of 20 wipers, 5 bags/case

# Sterile Pre-Wetted Wipers

STERILE



## One-Step Cleaning

A convenient and reliable alternative to wipers and alcohol spray bottles, Texwipe Sterile Pre-Wetted Wipers consist of cleanroom wipers with an optimum level of 70% isopropyl alcohol, gamma irradiated to a  $10^{-6}$  sterility assurance level.

- **Optimum saturation** — wet enough to clean surfaces and dry enough to remove contamination
- **Easier** — eliminates IPA storage, blending, filtration, in-house sterilization and validation documentation
- **Safer** — reduces spill potential, reduces solvent disposal, eliminates overspray of IPA and reduces VOCs

Sterile Pre-Wetted Wipers are ideal for cleaning aseptic fill areas, biotechnology production equipment, microbiological labs and isolators.

All Texwipe Sterile Pre-Wetted Wipers feature:

- **Gamma-irradiated to a  $10^{-6}$  sterility assurance level**
- **Clearly visible lot code and expiration date**
- **Lot-specific Certificates of Irradiation and Compliance**







### ISO Class 3-5



#### Sterile AlphaSat® 10 with WFI

SterileWipe™ LP 10 sealed-border, polyester wiper pre-wetted with 70% IPA/30% WFI. Ideal for cleaning equipment and environmental surfaces in critical aseptic environments.

Product Number	Description	Packaging
<b>TX3280</b>	12" x 12" nominal (31 cm x 31 cm)	50 wipers/flatpack, 5 flatpacks/case



#### Sterile AlphaSat® 10

SterileWipe™ LP 10 sealed-border, folded polyester wiper pre-wetted with 70% IPA. Packaged in a small, resealable flexpack sized appropriately for isolators and biological safety cabinets.

Product Number	Description	Packaging
<b>TX3285</b>	9" x 9" nominal (23 cm x 23 cm)	20 wipers/flexpack, 20 flexpacks/case

### ISO Class 5-6



#### Sterile AlphaSat® with AlphaSorb® HC

Sterile AlphaSorb® HC standard edge, polyester wiper pre-wetted with 70% IPA. Excellent for cleaning equipment and environmental surfaces in aseptic environments. Use for final wipedown to remove residual disinfectant.

Product Number	Description	Packaging
<b>TX3252</b>	12" x 12" nominal (31 cm x 31 cm)	25 wipers/flatpack, 5 flatpacks/case

### ISO Class 6-7



#### Sterile TechniSat®

SterileWipe™ HS II cellulose/polyester folded wiper pre-wetted with 70% IPA. All-purpose, durable pre-wetted wiper.

Product Number	Description	Packaging
<b>TX3214</b>	9" x 11" (23 cm x 28 cm)	50 wipers/flexpack, 24 flexpacks/case
<b>TX3217</b>	9" x 11" (23 cm x 28 cm)	20 wipers/flexpack, 24 flexpacks/case



#### Sterile PolySat®

Polypropylene folded wiper pre-wetted with 70% IPA. Economical, all-purpose pre-wetted wiper.

Product Number	Description	Packaging
<b>TX3213</b>	9" x 11" (23 cm x 28 cm)	50 wipers/flexpack, 24 flexpacks/case
<b>TX3216</b>	9" x 11" (23 cm x 28 cm)	20 wipers/flexpack, 24 flexpacks/case

## The First Name in Wipers



Texwipe cleanroom wipers efficiently clean, wipe and absorb spills in all classes of controlled environments. They also function as disposable work surfaces for bench top operations. The product family offers a range of properties and characteristics suitable for most applications.

### Polyester Knit Wipers

The polyester knit wipers offer the highest level of purity for ISO Class 3-7 environments. Their superior strength is ideal for rigorous wiping. They are compatible with detergents, disinfectants, alcohol, bleach and sterilants.

### Nonwoven Wipers

The nonwoven wipers consist of hydroentangled polyester and cellulose fibers. They are much cleaner and stronger than paper towels, and more economical than polyester knit wipers. They are ideal for all controlled environments ISO Class 7 and higher.

## ISO Class 1-3

## Polyester Knit Wipers



### AlphaSorb® 10

Premium two-ply, sealed-border polyester wiper for high sorption. Abrasion-resistant wiper with soft, flexible hand, designed to pick up large amounts of liquid and to apply chemicals. Ideal for use when one wiper is needed for multiple environments.

Product Number	Description	Packaging
TX1052	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 2 inner bags of 50, 6 bags/case



### TexTra™ 10

Sealed-border, heavy-weight polyester wiper for high sorption. Single-ply wiper with the liquid-holding capacity of a two-ply wiper. Good abrasion resistance, ideal for use when one wiper is needed for multiple environments.

Product Number	Description	Packaging
TX2452	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 2 inner bags of 50 wipers, 6 bags/case

## ISO Class 4-5

## Polyester Knit Wipers



### AlphaWipe®

Premium cut-edge, heavy-weight all-purpose polyester wiper for critical environments.

Product Number	Description	Packaging
<b>TX1003</b>	11" x 12" nominal (28 cm x 31 cm)	250 wipers/bag, 6 bags/case
<b>TX1013</b>	12" x 12" nominal (31 cm x 31 cm)	75 wipers/bag, double bagged, 10 bags/case



### TexTra™

Cut-edge heavy-weight, single-ply polyester wiper with sorption similar to a two-ply polyester wiper and good abrasion resistance.

Product Number	Description	Packaging
<b>TX2424</b>	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 2 inner bags of 50, 8 bags/case



### AlphaSorb® HC

Large, cut-edge two-ply polyester wiper for highest sorption capacity.

Product Number	Description	Packaging
<b>TX2412</b>	12" x 12" nominal (31 cm x 31 cm)	100 wipers/bag, 2 inner bags of 50 wipers, 6 bags/case
<b>TX2418</b>	18" x 18" nominal (46 cm x 46 cm)	100 wipers/bag, 2 inner bags of 50 wipers, 3 bags/case

## ISO Class 5, 6-9

## Nonwoven Wipers



### TechniCloth®

TechniCloth offers the lowest level of releasable ions and non-volatile contaminants; contains no chemical binders and tolerates most standard solvents. Its blend of 55% cellulose and 45% polyester provides the optimum balance of sorption, durability and cleanliness.

Product Number	Description	Packaging
<b>TX612</b>	12" x 12" nominal (31 cm x 31 cm)	150 wipers/bag, double bagged, 10 bags/case



### TechniCloth® II

TechniCloth II features the proprietary ULP (ultra-low particulate) treatment. It releases significantly fewer particles and fibers than the standard TechniCloth wiper.

Product Number	Description	Packaging
<b>TX1112</b>	12" x 12" nominal (31 cm x 31 cm)	150 wipers/bag, double bagged, 10 bags/case
<b>TX1118</b>	18" x 18" nominal (46 cm x 46 cm)	75 wipers/bag, double bagged, 10 bags/case



### VersaWipe®

Economical, lighter weight version of TechniCloth. It contains no chemical binders and tolerates most standard solvents.

Product Number	Description	Packaging
<b>TX622</b>	12" x 12" nominal (31 cm x 31 cm)	150 wipers/bag, double bagged, 10 bags/case


# Pre-Wetted Cleanroom Wipers

## Cleanroom Wipers with 70% IPA

Texwipe Pre-Wetted Wipers contain an optimum level of 70% isopropyl alcohol / 30% DIW solution. They are a convenient and reliable alternative to wipers and alcohol spray bottles:

- **Optimum saturation** — wet enough to clean surfaces and dry enough to remove contamination
- **Easier** — eliminates IPA storage, blending and filtration
- **Safer** — reduces spill potential, reduces solvent disposal; eliminates overspray of IPA, and reduces VOCs

The product family includes wipers designed for all classes of controlled environments. They are ideal for critical cleaning of pharmaceutical and biotechnology production equipment, microbiological labs and isolators.

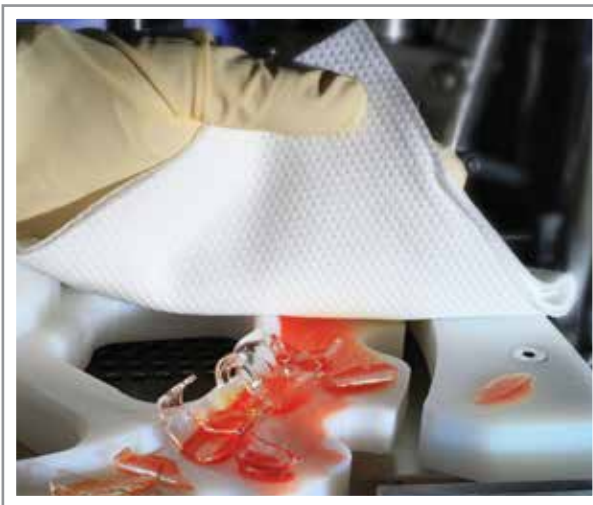

### ISO Class 1-3



#### QuanSat™ with Vectra® QuanTex™

Premium Polyester sealed border wiper for ultra-clean environments. Agarose Overlay and LAL testing reports available.

Product Number	Description	Packaging
<b>TX1084</b>	9" x 9" nominal (23 cm x 23 cm)	50 wipers/flatpack, 12 flatpacks/case

### ISO Class 4-5



#### AlphaSat® with AlphaWipe®

Heavy-weight polyester wiper for overall effective cleaning.

Product Number	Description	Packaging
<b>TX1034</b>	4" x 4" nominal (10 cm x 10 cm)	200 wipers/bucket, 4 buckets/case
<b>TX1036</b>	6" x 6" nominal (15 cm x 15 cm)	75 wipers/flip top tub, 12 fliptop tubs/case
<b>TX1039</b>	9" x 9" nominal (23 cm x 23 cm)	50 wipers/bucket, 4 buckets/case

### ISO Class 5



#### TechniSat®

Polyester/cellulose wiper for optimum cleaning and sorption.

Product Number	Description	Packaging
<b>TX1065</b>	9" x 11" nominal (23 cm x 28 cm)	50 wipers/flexpack, 24 flexpacks/case



#### PolySat®

Polypropylene wipers pre-wetted with 70% IPA

Product Number	Description	Packaging
<b>TX1051</b>	9" x 11" (22.9 cm x 27.9 cm)	50 wipers/flexpack, 24 flexpacks/case
<b>TX8723</b>	6" x 11" (15 cm x 28 cm)	75 wipers/ flexpack, 24 flexpacks/case

# Vertex® Pre-Wetted Wipers



### The Right Levels.

Overly saturated wipers can leave behind a pool of contaminants and insufficiently saturated wipers can perform poorly. Vertex Pre-Wetted Wipers contain just the right levels to clean contamination without over- or under-wetting surfaces.

### Same Wiper, Time After Time.

The Vertex guarantees that we produce highly consistent wipers time after time. No variation in particulates or fibers, just the highest level of assurance that you will receive a highly consistent wiper, time after time.

### Lowest Risk of Contamination.

From start to finish, this proprietary system isolates the wiper material from human contact. No people, no hands, no added risk of contamination.



### Vertex Pre-Wetted Wipers

Product Number	Description	Packaging
TX42P	Vertex 12" x 12" sealed-edge pre-wetted with 70% USP Purified IPA / 30% DIW wipers	50 wipers/bag, 4 bags/case
TX49P	Vertex 9" x 9" sealed-edge pre-wetted with 70% USP Purified IPA / 30% DIW wipers	75 wipers/bag, 4 bags/case

### STERILE



### Sterile Vertex Pre-Wetted Wipers

Product Number	Description	Packaging
TX3042P	Sterile Vertex 12" x 12" sealed-edge pre-wetted with 70% USP Purified IPA / 30% DIW wipers	25 wipers/bag, 5 bags/case
TX3044P	Sterile Vertex® 12" x 12" sealed-edge wipers pre-wetted with 70% USP-Grade Denatured Ethyl Alcohol and USP Purified Water	25 wipers/bag, 5 bags/case
TX3049P	Sterile Vertex 9" x 9" sealed-edge pre-wetted with 70% USP Purified IPA / 30% DIW wipers	25 wipers/bag, 5 bags/case





## Sterile Isopropyl Alcohol (IPA)

Sterile 70% Isopropanol contains 70% by volume USP-grade Isopropanol (isopropyl alcohol) and 30% USP-purified water.

Available in trigger spray and bulk containers. The trigger sprays are fully assembled with a ready-to-use applicator; no unpacking, handling or attachment of a separate spray head is required.

The solution is submicron filtered, filled into cleaned containers, double-bagged and gamma-irradiated. Each bottle is lot coded with an expiration date for easy recordkeeping.

### Sterile 70% Isopropanol

Specifically formulated alcohols for high-performance cleaning in sterile environments

Product Number	Description	Packaging
TX8270	Trigger spray, 8 fluid ounces (237 ml)	12 polybottles/case
TX3270	Trigger spray, 16 fluid ounces (473 ml)	12 polybottles/case
TX3273	Trigger spray, 32 fluid ounces (946 mL)	12 polybottles/case
TX3290	Bulk container, 1 gallon (3.8 liters).	4 polybottles/case

# Quaternary Ammonium Disinfectant



# 25% More Kill Claims



TX650  
Ready-to-Use Trigger Spray

## TexQ™ Disinfectant

### One-Step Cleaner and Disinfectant

TexQ Disinfectant is an EPA-registered quaternary ammonium compound (quat) which serves as an effective one-step cleaner/disinfectant. It is effective against a broad spectrum of bacteria, virus and fungi, and inhibits the growth of mold and mildew and their odors when used as directed.

**1 Qualification, 2 Solutions.** TexQ is available in trigger spray and concentrate formulas - both with the same active ingredients and final concentration. Qualify the product once and use either form with confidence as the end-use dilution at point-of-use is the same for both solutions.

#### TexQ Disinfectant

Product Number	Description	Packaging
TX650	TexQ Disinfectant Ready-to-Use (RTU)	22 oz. trigger spray, 12 bottles/case
TX651	TexQ Disinfectant Concentrate	1 gallon concentrate, 4 bottles/case

# Cleanroom Mops

## From Floor to Ceiling

Texwipe Cleanroom Mops combine the effectiveness of cleanroom wipers with convenient handles to wipe large surfaces and difficult-to-reach areas.

Our products include large flat mops for wiping large environmental surfaces, small flat mops for hard-to-reach areas, and large volume mops for floors.

Also available are mop covers and replacement parts for the Texwipe Cleanroom Mop products.



TX7106  
BetaMop™ II



### BetaMop™ II

15-ounce tubular polyester mop for efficient and economical floor cleaning in critical or clean environments. All components are compatible with autoclaves, chemicals, VHP or EO.

Product Number	Name / Description	Packaging
TX7106	BetaMop™ II Mop with head refill	1 mop/case

### Mop Covers & Refills

TX716R	BetaMop™ II Refill heads for use with TX106 BetaMop II	6 heads/case
--------	---	--------------

TX7102  
ClipperMop™



### ClipperMop™

Mop systems that use 9" x 9" or 12" x 12" wipers. The entire system tolerates sterilization by chemicals, VHP or EO.

Product Number	Description	Packaging
TX7102	ClipperMop™ 7" mop head, extension handle	1 mop/case
TX7103	ClipperMop™ 11" mop head, extension handle	1 mop/case

### Mop Accessories & Replacement Parts

TX7102C	ClipperMop™ Clips Polyethylene tether clips and nylon screws for use with TX7102 and TX7103	4 clips/bag = 80 clips/cs 2 screws/bag = 40 screws/cs 20 bags/case
TX7111A	ClipperMop™ Replacement Pads (7102) Foam pads and white fasteners for use with TX7102	10 pads/bag = 100 pads/cs 4 fasteners/bag = 40 fasteners 10 bags/case
TX7111B	ClipperMop™ Replacement Pads (7103) Foam pads and white fasteners for use with TX7103	10 pads/bag = 100 pads/cs 4 fasteners/bag = 40 fasteners 10 bags/case



TX7108  
AlphaMop™



TX7118F  
AlphaMop™ Foam Pad Replacement Kit

### AlphaMop™

Ideal for applications requiring low particle, extractable and ionic contamination. Superior for wall, floor and ceiling cleaning in ultraclean manufacturing environments. Autoclavable.

Product Number	Name / Description	Packaging
<b>TX7108</b>	AlphaMop™ Cleanroom mop with handle and single head assembly	1 mop/case
<b>TX7108A</b>	Autoclavable AlphaMop™ 15" x 8" thermoplastic single head assembly 60" fiberglass handle	1 mop/case

#### Mop Covers & Refills

<b>TX7026</b>	AlphaSat® Pre-Wetted Mop Covers 8" x 15" (20 cm x 38 cm) pre-wetted with 6% IPA	25 mop covers/bag 4 bags/case = 100 mop covers 4 foam pads/case
<b>TX7118</b>	AlphaMop™ Mop Covers 8" x 15" (20 cm x 38 cm) Polyester mop covers and foam pads	150 mop covers/case (6 bags of 25 mop covers) 6 foam pads/case
<b>TX7118F</b>	AlphaMop™ Foam Pad Replacement Kit Foam pads and clips for use with AlphaMop™ TX7108 and TX7108A	25 foam pads and 4 clips/bag 3 bags/case
<b>TX3268</b>	Sterile AlphaMop™ Covers polyester replacement covers and pads	125 covers and 25 pads/case (25 bags of 5 covers and 1 pad)

#### Mop Accessories & Replacement Parts

<b>TX7160</b>	Autoclavable Extension Handle 60" - 112"	1 handle/case
---------------	---	---------------

### Mini AlphaMop™

Mop system for cleaning isolators, biosafety cabinets, glove boxes and laminar flow hoods. Compatible with sterilization by autoclave, chemicals, VHP or EO.

Product Number	Name / Description	Packaging
<b>TX7101</b>	Isolator Cleaning Tool with Mini AlphaMop™ 3.7" x 7.1" (9.4 cm x 18.1 cm) Mini cleanroom mop with 2 extension handles plus foam pad and starter pack of covers	1 mop/case 2 handles (12" and 18") 1 foam pad 6 covers (TX7114)
<b>TX7104</b>	Mini AlphaMop™ 3.7" x 7.1" (9.4 cm x 18.1 cm) Mini cleanroom mop with 1 telescoping handle plus foam pad and starter pack of covers	1 mop/case 1 telescoping handle (telescopes from 2' 5" to 4' 5") 1 foam pad 6 covers (TX7114)

#### Mop Covers & Refills

<b>TX3264</b>	Sterile Mini AlphaMop™ Covers polyester replacement covers and pads	120 covers and 12 pads/case (12 bags of 10 covers and 1 pad)
<b>TX7114M</b>	Mini AlphaMop™ Replacement Covers Microdenier polyester mop covers	150 mop covers/case (6 bags of 25 mop covers) 6 foam pads/case

#### Mop Accessories & Replacement Parts

<b>TX7160</b>	Autoclavable Extension Handle 60" - 112"	1 handle/case
---------------	---	---------------



TX7101  
Isolator Cleaning Tool with  
Mini AlphaMop™

# Cleaning Validation Swabs

## Green Handles are the Gold Standard

Texwipe CleanTips® Swabs, with their trademarked green handles, are the global standard for cleaning validation in FDA regulated industries.

CleanTips® Swabs are engineered for the stringent requirements of both specific (HPLC, IMS, UV-vis) and non-specific (TOC) analytical methods.

The HPLC and IMS swabs are solvent purified to reduce background interference for more accurate chemical analysis. These swabs are effective with a wide range of

diluents. The TOC sampling swabs are certified to <50 µg/L and are designed for the efficient collection and recovery of residual organic contamination.



### Cleaning Validation Swabs



#### Low TOC Sampling Swab (SnapSwab™)

Knitted polyester swab with broad paddle head purified and certified to <50 µg/L TOC. Snap off head feature. Also available in TOC Cleaning Validation Kits.

Product Number	Description	Packaging
<b>TX714K</b>	5.0" overall swab length	20 swabs/bag, 50 bags/case



#### Sampling Swab

Solvent cleaned knitted polyester swab optimized for HPLC, IMS or UV-vis analysis. Snap off head feature.

Product Number	Description	Packaging
<b>TX715</b>	5.0" overall swab length	100 swabs/bag (2 inner bags of 50 swabs), 10 bags/case



#### Large Sampling Swab

Hydroentangled polyester head significantly reduces background interference in HPLC, IMS or UV-vis analysis. Snap off head feature.

Product Number	Description	Packaging
<b>TX716</b>	5.0" overall swab length	20 swabs/bag, 50 bags/case



#### Low TOC Swab (with Long Handle)

Knitted polyester swab with long handle purified and certified to <50 µg/L TOC. Flexible to provide contact in difficult to reach sampling areas.

Product Number	Description	Packaging
<b>TX761K</b>	6.4" overall swab length	20 swabs/bag, 50 bags/case





## The Preferred Tools

Texwipe's CleanTips® Swabs are the preferred tools for rigorous cleaning in tight areas, particularly ideal for pharmaceutical and biotech work environments.

The Alpha® Series Swabs feature cleanroom-laundered knit-polyester heads for extraordinary cleanliness. The polyester heads are thermally molded onto the polypropylene handles, eliminating the potential of contamination from adhesives. The perimeter edges are sealed to prevent fiber release.

Alpha® Swabs can be sterilized by autoclave, EO or steam. They are compatible with a wide range of detergents, disinfectants and cleaning agents.



### Alpha® Swabs



#### Large Alpha® Swab

The large paddle-style head and sturdy handle for removing contaminants in a broad path.

Product Number	Description	Packaging
<b>TX714A</b>	5.0" overall swab length	100 swabs/bag (2 inner bags of 50 swabs), 10 bags/case



#### Small and Mini Alpha® Swabs

These short-handled, small-tipped swabs are designed for cleaning tight areas where rigidity on the sides of the head is required.



Product Number	Description	Packaging
<b>TX743B</b>	<b>Small Alpha Swab</b> 2.8" overall swab length	500 swabs/bag (5 inner bags of 100 swabs), 5 bags/case



Product Number	Description	Packaging
<b>TX754B</b>	<b>Mini Alpha Swab</b> 2.7" overall swab length	500 swabs/bag (5 inner bags of 100 swabs), 5 bags/case

#### Micro Alpha® Swab

The flexible tip enables cleaning at an angle perpendicular to the swab length in narrow areas.

Product Number	Description	Packaging
<b>TX758B</b>	2.8" overall swab length	500 swabs/bag (5 inner bags of 100 swabs), 5 bags/case



#### Alpha® Swab with Long Handle

The medium-size paddle-shaped tip and long, flexible handle make this swab the perfect tool for use in hard-to-reach and partially obscured areas.

Product Number	Description	Packaging
<b>TX761</b>	6.4" overall swab length	100 swabs/bag (1 inner bag of 100 swabs), 10 bags/case



*Additional information including dimensions for these and more than 30 other CleanTips Swabs is available at [www.texwipe.com](http://www.texwipe.com).*

# TOC Cleaning Validation Kits



TX3342  
Bulk TOC Cleaning Validation Kit



## TOC Cleaning Validation Kits

Texwipe created TOC (Total Organic Carbon) Cleaning Validation Kits to simplify cleaning validation by TOC analysis. The kits include swabs, vials and labels in a cleanroom-compatible polypropylene box.

The vials are certified to <10 µg/L and have a 40 mL capacity. The sampling swabs are certified to <50 µg/L and

are designed for the efficient collection and recovery of residual organic contamination. The TX714K SnapSwab™ has notches for easy breaking.

The polypropylene box is compatible with cleanrooms; it protects the integrity of the vials both before and after sampling; and it offers safe transportation of the glass vials.

### TOC Cleaning Validation Kits

Product Number	Description	Packaging
TX3340	TOC Cleaning Validation Kit, 12 vials + 24 swabs + 2 labels	18 kits/case
TX3342	Bulk TOC Cleaning Validation Kit, 72 vials + 144 swabs + 72 labels	1 kit/case

## Leave your mark

Texwipe's quick-ship pressure sensitive cleanroom labels are designed for critical environments and are made from material that is low out-gassing, has ultra-low leachables and low particulates. Each label is tested for surface particles, cleaned, double bagged and packaged in a controlled environment.



### Blank Cleanroom Labels

Product #	Stock Type	Size	Description	Packaging
<b>Autoclave</b>				
<b>TX515A</b>	Autoclave	1.5" round	Autoclave Labels, Blank, 1.5" round - Autoclave adhesive	500 labels/roll, 4 rolls/case
<b>Direct Thermal</b>				
<b>TX515D</b>	Direct Thermal	1.5" x 1.5"	Direct Thermal Labels, Blank, 1.5" x 1.5" - Permanent adhesive	1000 labels/roll, 4 rolls/case
<b>TX532D</b>	Direct Thermal	3" x 2"	Direct Thermal Labels, Blank, 3" x 2" - Permanent adhesive	1000 labels/roll, 4 rolls/case
<b>TX542D</b>	Direct Thermal	4" x 2"	Direct Thermal Labels, Blank, 4" x 2" - Permanent adhesive	1000 labels/roll, 4 rolls/case
<b>Gloss</b>				
<b>TX515G</b>	Gloss	1.5" x 1.5"	Gloss Labels, Blank, 1.5" x 1.5" - Removable adhesive	1000 labels/roll, 4 rolls/case
<b>TX532G</b>	Gloss	3" x 2"	Gloss Labels, Blank, 3" x 2" - Removable adhesive	1000 labels/roll, 4 rolls/case
<b>TX542G</b>	Gloss	4" x 2"	Gloss Labels, Blank, 4" x 2" - Permanent adhesive	1000 labels/roll, 4 rolls/case
<b>Matte</b>				
<b>TX515M</b>	Matte	1.5" x 1.5"	Matte Labels, Blank, 1.5" x 1.5" - Removable adhesive	1000 labels/roll, 4 rolls/case
<b>TX532M</b>	Matte	3" x 2"	Matte Labels, Blank, 3" x 2" - Removable adhesive	1000 labels/roll, 4 rolls/case
<b>TX542M</b>	Matte	4" x 2"	Matte Labels, Blank, 4" x 2" - Removable adhesive	1000 labels/roll, 4 rolls/case

### Custom Cleanroom Labels

Texwipe specializes in printing custom cleanroom labels. Please contact your Regional Sales Manager for more information.

### Pre-Printed Cleanroom Labels

Product Number	Name / Description	Packaging
<b>TX522A</b> (formerly ALAU275)	Pre-Printed Sterilization Indicator, Autoclave 2" x 2" - Autoclave adhesive	500 labels/roll 1 roll/case

## Documented Clean

Releasing significantly fewer particles and fibers during copying, printing and handling than conventional bond paper. TexWrite® paper and notebooks are engineered for critical environments.

### Loose-Leaf Sheets.

Available in **premium synthetic** (TX5800) and **saturated cellulose** papers in 30#, 22# and 18# weights. Latex-saturated sheets are available in blue, green, yellow and white.

Premium synthetic sheets are available in a pleasing blue-green color and are compatible with sterilization by gamma irradiation, EO or steam.

### Journals & Notebooks.

Ideal for documentation and storage of notes and records in ISO Class 4 - 6 cleanrooms. Our proprietary Bag-Within-A-Bag® technology helps guard against contamination by allowing alcohol wipe down (of exterior bag) before introduction of interior bag into cleanrooms, sterile suites and aseptic environments.

### Binders & Clipboards.

Designed for cleanroom equipment manuals, operating procedures and work instructions, the TexWrite® **three ring binder** features a durable polyethylene cover with a plastic ring mechanism and the **Cleanroom Clipboard** is the ideal portable writing surface.



TX5820  
TexWrite®  
Adhesive Pads



TX5835  
TexWrite® Clipboard



TX5708  
Cleanroom  
Notebook



Cleanroom  
Spiral Notebooks

### Adhesive Note Pads

Product Number	Name / Description	Packaging
TX5820	TexWrite® Adhesive Pads 3" x 4" (76 mm x 102 mm) 22 lb. blue paper stock, unlined	72 sheets/pad 10 pads/box (case) 720 sheets

### Binders & Clipboards


Product Number	Name / Description	Packaging
TX5810	Three-Ring Binder, TexWrite® 1" (25 mm) capacity slide ring mechanism	10 binders/case
TX5835	TexWrite® Clipboard 9" x 13"	12 clipboards/case

### Journals & Notebooks

Product Number	Name / Description	Packaging
TX5708	Cleanroom Notebook, TexWrite® 8½" x 11" (216 mm x 279 mm) Blue paper stock, college ruled, numbered pages	10 notebooks/box (case)
TX5740	Cleanroom Spiral Notebook, TexWrite® 22 8½" x 11" (216 mm x 279 mm) 22 lb. white paper stock, college ruled	10 notebooks/box (case)
TX5741	Cleanroom Spiral Notebook, TexWrite® 22 5½" x 8½" (140 mm x 216 mm) 22 lb. white paper stock, college ruled	10 notebooks/box (case)
TX5742	Cleanroom Spiral Notebook, TexWrite® 22 3" x 5" (76 mm x 127 mm) 22 lb. white paper stock, college ruled	20 notebooks/box (case)



**TX5800**  
TexWrite® MP 10  
synthetic sheet

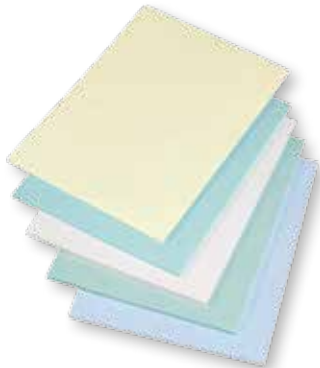

TexWrite® Paper

## Loose Leaf Sheets

Product Number	Name / Description	Packaging
<b>TX5800</b>	TexWrite® MP 10 8½" x 11" (216 mm x 279 mm) white/blue-green 100% synthetic sheet, unlined	100 sheets/pack 5 packs/box = 500 sheets
<b>TX5812</b>	TexWrite® 22 Loose Sheets, BLUE 8½" x 11" (216 mm x 279 mm) 22 lb. blue paper stock, unlined	250 sheets/pack 10 packs/box (case) 2,500 sheets
<b>TX5814</b>	TexWrite® 22 Loose Sheets, YELLOW 8½" x 11" (216 mm x 279 mm) 22 lb. yellow paper stock, unlined	250 sheets/pack 10 packs/box (case) 2,500 sheets
<b>TX5815</b>	TexWrite® 22 Loose Sheets, WHITE 8½" x 11" (216 mm x 279 mm) 22 lb. white paper stock, unlined	250 sheets/pack 10 packs/box (case) 2,500 sheets
<b>TX5816</b>	TexWrite® 22 3-Hole Punch Sheets, BLUE 8½" x 11" (216 mm x 279 mm) 22 lb. blue paper stock, 3-hole punched, unlined	250 sheets/pack 10 packs/box (case) 2,500 sheets
<b>TX5831</b>	TexWrite® 22 Loose Sheets, GREEN 8½" x 11" (216 mm x 279 mm) 22 lb. green paper stock, unlined	250 sheets/pack 10 packs/box (case) 2,500 sheets
<b>TX5832</b>	TexWrite® 30 Loose Sheets, WHITE 8½" x 11" (216 mm x 279 mm) 30 lb. white paper stock, unlined	250 sheets/pack 7 packs/box (case) 1,750 sheets
<b>TX5862</b>	TexWrite® 18 Loose Sheets, BLUE 8½" x 11" (216 mm x 279 mm) 18 lb. blue paper stock, unlined	250 sheets/pack 10 packs/box (case) 2,500 sheets
<b>TX5916</b>	TexWrite® 22 3-Hole Punch Sheets, WHITE 8½" x 11" (216 mm x 279 mm) 22 lb. white paper stock, 3-hole punched, unlined	250 sheets/pack 10 packs/box (case) 2,500 sheets



# Cleanroom Tapes

## Cleanroom Tapes

Texwipe offers a wide range of Cleanroom Tapes, constructed from the highest quality materials. Available in a variety of sizes, colors and materials, we have a solution to your specific application.

Full line of tape constructions suitable for most applications and environments. We selected specifically engineered tape choices to suit an array of preferences, protocols and needs.






















**LDPE with Acrylic** – strong adhesion, hand-tearable, halogen and leachable construction

**LDPE with Synthetic Rubber** – hand-tearable, low halogens and leachables





















**PVC with Synthetic Rubber** – extremely conformable tape with aggressive adhesion, with clean peel from most cleanroom surfaces























### LDPE with Acrylic Adhesive (TPA Series)

Product Number	Color	Size	Packaging
TPA0596CL		0.5" x 36 yds	24 bags/4 rolls ea. / 96 rolls total
TPA0596WH		0.5" x 36 yds	24 bags/4 rolls ea. / 96 rolls total
TPA0596BL		0.5" x 36 yds	24 bags/4 rolls ea. / 96 rolls total
TPA0596GR		0.5" x 36 yds	24 bags/4 rolls ea. / 96 rolls total
TPA0596OR		0.5" x 36 yds	24 bags/4 rolls ea. / 96 rolls total
TPA1048CL		1" x 36 yds	24 bags/2 rolls ea. / 48 rolls total
TPA1048WH		1" x 36 yds	24 bags/2 rolls ea. / 48 rolls total
TPA1048BL		1" x 36 yds	24 bags/2 rolls ea. / 48 rolls total
TPA1048GR		1" x 36 yds	24 bags/2 rolls ea. / 48 rolls total
TPA1048OR		1" x 36 yds	24 bags/2 rolls ea. / 48 rolls total
TPA2024CL		2" x 36 yds	24 rolls/case
TPA2024WH		2" x 36 yds	24 rolls/case
TPA2024BL		2" x 36 yds	24 rolls/case
TPA2024GR		2" x 36 yds	24 rolls/case
TPA2024OR		2" x 36 yds	24 rolls/case
TPA3016CL		3" x 36 yds	16 rolls/case
TPA3016WH		3" x 36 yds	16 rolls/case
TPA3016BL		3" x 36 yds	16 rolls/case
TPA3016GR		3" x 36 yds	16 rolls/case
TPA3016OR		3" x 36 yds	16 rolls/case
TPA2036WHC		Cuff Sealer Cleanroom Tape White, Perforated / 2" wide	24 rolls/case


## LDPE with Synthetic Rubber Adhesive (TPR Series)

Product Number	Color	Size	Packaging
TPR0596YE		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TPR0596WH		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TPR0596BL		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TPR0596GR		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TPR0596RD		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TPR1048YE		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TPR1048WH		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TPR1048BL		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TPR1048GR		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TPR1048RD		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TPR2024YE		2" x 36 yds	24 rolls/case
TPR2024WH		2" x 36 yds	24 rolls/case
TPR2024BL		2" x 36 yds	24 rolls/case
TPR2024GR		2" x 36 yds	24 rolls/case
TPR2024RD		2" x 36 yds	24 rolls/case
TPR3016YE		3" x 36 yds	16 rolls/case
TPR3016WH		3" x 36 yds	16 rolls/case
TPR3016BL		3" x 36 yds	16 rolls/case
TPR3016GR		3" x 36 yds	16 rolls/case
TPR3016RD		3" x 36 yds	16 rolls/case

## PVC with Synthetic Rubber Adhesive (TVR Series)

Product Number	Color	Size	Packaging
TVR0596YE		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TVR0596WH		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TVR0596BL		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TVR0596GR		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TVR0596RD		0.5" x 36 yds	24 bags/4 rolls ea. - 96 rolls total
TVR1048YE		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TVR1048WH		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TVR1048BL		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TVR1048GR		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TVR1048RD		1" x 36 yds	24 bags/2 rolls ea. - 48 rolls total
TVR2024YE		2" x 36 yds	24 rolls/case
TVR2024WH		2" x 36 yds	24 rolls/case
TVR2024BL		2" x 36 yds	24 rolls/case
TVR2024GR		2" x 36 yds	24 rolls/case
TVR2024RD		2" x 36 yds	24 rolls/case
TVR3016YE		3" x 36 yds	16 rolls/case
TVR3016WH		3" x 36 yds	16 rolls/case
TVR3016BL		3" x 36 yds	16 rolls/case
TVR3016GR		3" x 36 yds	16 rolls/case
TVR3016RD		3" x 36 yds	16 rolls/case

# Cleanroom Mats

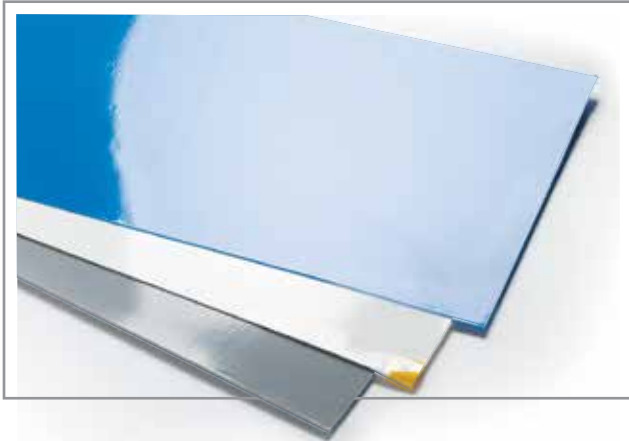

## Time-Saving, Floor-Saving

The CleanStep™ Adhesive Mat removes traffic-borne particulates prior to entering the clean work area. It features an antimicrobial additive in an environmentally friendly adhesive.

Whether cleaning protocols require portability or the surface just isn't conducive to adhesives, CleanStep's mix of sizes and colors are sure to fit the bill. Offering 36 mat/frame combinations, the easy selection will improve your environment and save you time.

### 18 in. x 36 in.

Available in 30 and 60 layer and in Blue, Gray, and White

Product Number	Description	Packaging
<b>AMA183681B</b>	Blue base, clear film / 18" x 36" / 30 layers	8 mats per case
<b>AMA183681G</b>	Gray base, clear film / 18" x 36" / 30 layers	8 mats per case
<b>AMA183681W</b>	White base, clear film / 18 in. x 36" / 30 layers	8 mats per case
<b>AMA183682B</b>	Blue Base, Clear Film / 18" x 36" / 60 layers	4 mats per case
<b>AMA183682G</b>	Gray Base, Clear Film / 18" x 36" / 60 layers	4 mats per case
<b>AMA183682W</b>	White Base, Clear Film / 18" x 36" / 60 layers	4 mats per case

### 25 in. x 45 in.

Available in 30 and 60 layer and in Blue, Gray, and White.

Product Number	Description	Packaging
<b>AMA254581B</b>	Blue base, clear film / 25" x 45" / 30 layers	8 mats per case
<b>AMA254581G</b>	Gray base, clear film / 25" x 45" / 30 layers	8 mats per case
<b>AMA254581W</b>	White base, clear film / 25" x 45" / 30 layers	8 mats per case
<b>AMA254582B</b>	Blue Base, Clear Film / 25" x 45" / 60 layers	4 mats per case
<b>AMA254582G</b>	Gray Base, Clear Film / 25" x 45" / 60 layers	4 mats per case
<b>AMA254582W</b>	White Base, Clear Film / 25" x 45" / 60 layers	4 mats per case

### 36 in. x 46 in.


Available in 30 and 60 layer and in Blue, Gray, and White.

Product Number	Description	Packaging
<b>AMA364681B</b>	Blue base, clear film / 36" x 46" / 30 layers	8 mats per case
<b>AMA364681G</b>	Gray base, clear film / 36" x 46" / 30 layers	8 mats per case
<b>AMA364681W</b>	White base, clear film / 36" x 46" / 30 layers	8 mats per case
<b>AMA364682B</b>	Blue Base, Clear Film / 36" x 46" / 60 layers	4 mats per case
<b>AMA364682G</b>	Gray Base, Clear Film / 36" x 46" / 60 layers	4 mats per case
<b>AMA364682W</b>	White Base, Clear Film / 36" x 46" / 60 layers	4 mats per case

### Frames

Frame replacement parts for CleanStep Mat system.

Product Number	Description	Packaging
<b>FB1937</b>	Frame: 19" x 37" (for mat size 18" x 36") Color: Black	1 frame per case
<b>FW1937</b>	Frame: 19" x 37" (for mat size 18" x 36") Color: White	1 frame per case
<b>FB2646</b>	Frame: 26" x 46" (for mat size 25" x 45") Color: Black	1 frame per case
<b>FW2646</b>	Frame: 26" x 46" (for mat size 25" x 45") Color: White	1 frame per case
<b>FB3747</b>	Frame: 37" x 47" (for mat size 36" x 46") Color: Black	1 frame per case
<b>FW3747</b>	Frame: 37" x 47" (for mat size 36" x 46") Color: White	1 frame per case

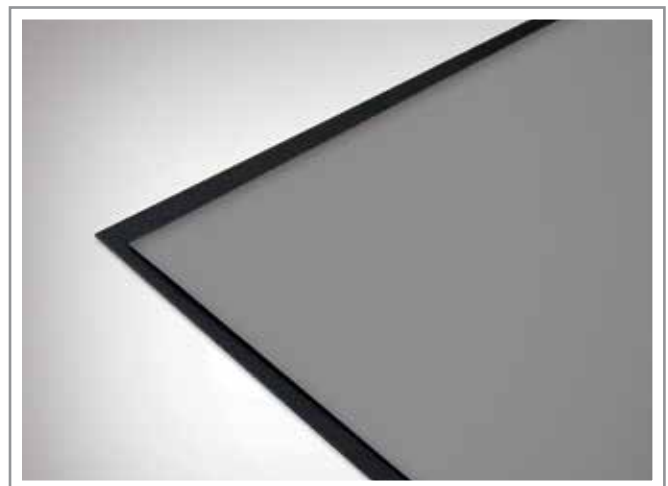

## Product Index

<b>Adhesive Mats</b>		Page
CleanStep™ Adhesive Mats and Frames		18

<b>Adhesive Note Pads</b>		Page
TX5280		20

<b>Cleanroom Mops</b>		Page
TX716R	BetaMop™ II Refills	17
TX7101	AlphaMop™ Isolator Cleaning Tool	17
TX7102	ClipperMop™	17
TX7103	ClipperMop™	17
TX7106	BetaMop™ II	17
TX7111	TechniCloth® Replacement Covers	17
TX7114	Polyester Replacement Covers	17

<b>Cleanroom Paper</b>		Page
TX5707	TexWrite® 22 Permanent-Bound Notebook	20
TX5708	TexWrite® 22 Permanent-Bound Notebook	20
TX5740	TexWrite® 22 Spiral-Bound Notebook	20
TX5741	TexWrite® 22 Spiral-Bound Notebook	20
TX5742	TexWrite® 22 Spiral-Bound Notebook	20
TX5800	Loose-Leaf Sheets, Synthetic	20
TX5812	Loose-Leaf Sheets, Latex/Cellulose	20
TX5815	Loose-Leaf Sheets, Latex/Cellulose	20
TX5835	TexWrite® Clipboard	20

<b>Cleanroom Swabs</b>		Page
TX714A	Large Alpha® Swab	16
TX743B	Small Alpha® Swab	16
TX754B	Mini Alpha® Swab	16
TX758B	Micro Alpha® Swab	16
TX761	Alpha® Swab with Long Handle	16

<b>Cleanroom Swabs - Cleaning Validation</b>		Page
TX714K	Alpha® Low TOC Sampling Swab	14
TX715	Alpha® Sampling Swab	14
TX716	Absorbond® Sampling Swab	14
TX761K	Low TOC Alpha® Swab with Long Handle	14

<b>Cleanroom Tapes</b>		Page
KÓGÓ Cleanroom Tapes		19

<b>Disinfectants</b>		Page
TX650	TexQ™ Disinfectant, Ready-to-Use™	13
TX651	TexQ™ Disinfectant, Concentrate™	13

### STERILE

<b>Sterile IPA</b>		Page
TX3270	Trigger Spray, Sterile IPA	12
TX3273	Trigger Spray, Sterile IPA	12
TX3290	Bulk Container, Sterile IPA	12
TX8270	Trigger Spray, Sterile IPA	12

<b>TOC Cleaning Validation Kits</b>		Page
TX3340	TOC Cleaning Validation Kit	15
TX3342	TOC Cleaning Validation Bulk Kit	15

## WIPERS

<b>Dry Wipers</b>		Page
TX22	Vertex® High Durability	2
TX29	Vertex® High Durability	2
TX42	Vertex® High Sorption	2
TX49	Vertex® High Sorption	2
TX59	Vertex® Microdenier	2
TX612	TechniCloth®	9
TX622	VersaWipe®	9
TX1013	AlphaWipe®	9
TX1052	AlphaSorb® 10	8
TX1112	TechniCloth® II	9
TX1118	TechniCloth® II	9
TX2412	AlphaSorb® HC	9
TX2418	AlphaSorb® HC	9
TX2424	TexTra™	9
TX2452	TexTra™ 10	8

### Pre-Wetted Wipers

TX42P	Vertex® Pre-Wetted	9
TX49P	Vertex® Pre-Wetted	9
TX1034	AlphaSat® with AlphaWipe®	8
TX1036	AlphaSat® with AlphaWipe®	8
TX1039	AlphaSat® with AlphaWipe®	8
TX1051	PolySat®	8
TX1065	TechniSat®	8
TX1084	QuanSat™ with Vectra® QuanTex™	8

### STERILE

<b>Sterile Dry Wipers</b>		Page
TX3042	Vertex® Sterile	3
TX3049	Vertex® Sterile	3
TX3052	Vertex® Sterile	3
TX3059	Vertex® Sterile	3
TX3210	SterileWipe™ HS II	5
TX3211	SterileWipe™ LP	5
TX3212	SterileWipe™ LP 10	5
TX3215	SterileWipe™ AS 10	5
TX3220	Sterile TexTra™	5
TX3224	Sterile TexTra™ 10	5
TX3225	Sterile TexTra™ 10	5

<b>Sterile Pre-Wetted Wipers</b>		Page
TX3042P	Vertex® Sterile Pre-Wetted	11
TX3044P	Vertex® Sterile Pre-Wetted	11
TX3049P	Vertex® Sterile Pre-Wetted	11
TX3213	Sterile PolySat®	7
TX3214	Sterile TechniSat®	7
TX3216	Sterile PolySat®	7
TX3217	Sterile TechniSat®	7
TX3252	Sterile AlphaSat® with AlphaSorb® HC	7
TX3280	Sterile AlphaSat® 10 with WFI	7
TX3285	Sterile AlphaSat® 10	7

*For more information contact us at the numbers below or [www.texwipe.com](http://www.texwipe.com)*

**North America Office** +1.336.996.7046

**Europe/Middle East Office** +31.88.1307.410

**Asia/Pacific Office** +65 6468 9433



Distributed in Australia & New Zealand by:

**ThermoFisher**  
SCIENTIFIC

**Australia:** For customer service, call 1300-735-292

Visit us online at: [www.thermofisher.com.au](http://www.thermofisher.com.au)

**New Zealand:** For customer service, call 0800-933-966

Visit us online at: [www.thermofisher.co.nz](http://www.thermofisher.co.nz)

©2015 Thermo Fisher Scientific Inc. All rights reserved.