

Nunculon Supra Surface for primary cells

Cell culture



Nunc cultureware can help you culture with confidence

Nunculon Supra Surface enables easier seeding and handling and improved confluency and cell yields over typical cell culture surfaces


Thermo Scientific™ Nunc™ cultureware products with the Thermo Scientific™ Nunculon™ Supra Surface feature a one-size-fits-all surface that has proven compatibility with >40 cell lines and cell types.

The Nunculon Supra Surface is designed to increase ease of use for many types of finicky cell lines and primary cells. It can also help:

- **Reduce time for adherence**, allowing for cells to start proliferating sooner
- **Improve morphology**, indicating healthier growing cells
- **Improve cell yield and confluency**, helping provide more cells in the same amount of culture time

Primary fibroblasts

Tissue culture surface

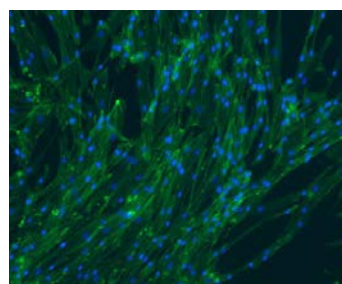

Nunculon Supra Surface


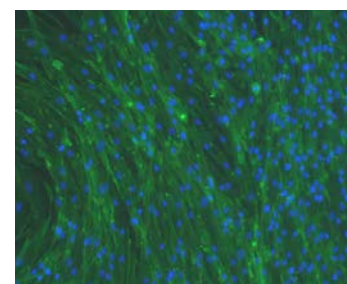


Figure 1. Adult human fibroblasts (HDFa) on tissue culture and Nunculon Supra surfaces. HDFa after immunofluorescence staining; f-actin is stained green and DNA is stained blue.


Primary endothelial cells

A

Tissue culture surface

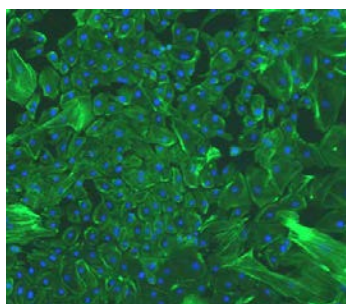


Nunculon Supra Surface

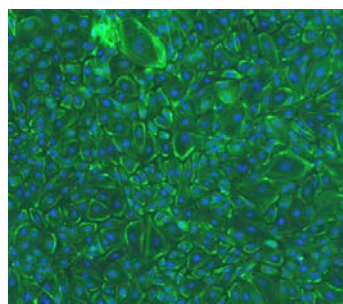

B

Tissue culture surface

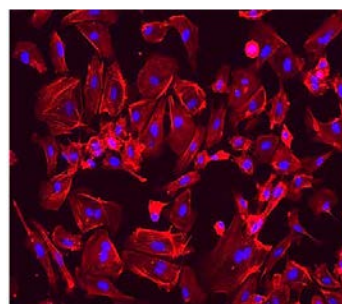

Nunculon Supra Surface

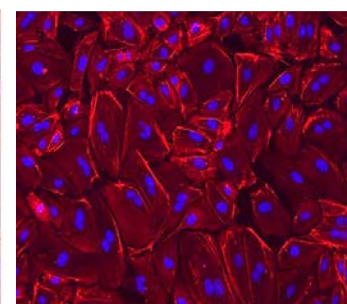


Figure 2. Two endothelial cell types cultured on tissue culture and Nunculon Supra surfaces. (A) Human aortic endothelial cells with immunofluorescence staining for f-actin (green) and DNA (blue). (B) Human umbilical vein endothelial cells with immunofluorescence staining for f-actin (red) and DNA (blue).

Ordering information

Type	Product	Cat. No.
Plates	Nunc Multidish with Nunclon Supra Surface, 6-well	140680
	Nunc Multidish with Nunclon Supra Surface, 12-well	140681
	Nunc Multidish with Nunclon Supra Surface, 24-well	140682
	Nunc Multidish with Nunclon Supra Surface, 48-well	140683
Dishes	Nunc EasYDish Dish with Nunclon Supra Surface, 35 mm	150470
	Nunc EasYDish Dish with Nunclon Supra Surface, 60 mm	150472
	Nunc EasYDish Dish with Nunclon Supra Surface, 100 x 15 mm	150474
	Nunc EasYDish Dish with Nunclon Supra Surface, 100 x 20 mm	150476
	Nunc EasYDish Dish with Nunclon Supra Surface, 150 mm	150478
Flasks	Nunc EasYFlask Flask with Nunclon Supra Surface, T25, filter cap	156372
	Nunc EasYFlask Flask with Nunclon Supra Surface, T25, solid cap	156374
	Nunc EasYFlask Flask with Nunclon Supra Surface, T75, filter cap	156376
	Nunc EasYFlask Flask with Nunclon Supra Surface, T75, solid cap	156378
	Nunc EasYFlask Flask with Nunclon Supra Surface, T175, filter cap	156380
	Nunc EasYFlask Flask with Nunclon Supra Surface, T175, solid cap	156382
96-well plates	Nunc Microwell 96-Well Microplate with Nunclon Supra Surface, with lid	167013
Cells and cell lines	Gibco Human Dermal Fibroblasts, adult (HDFa)	C0135C
	Gibco Human Aortic Endothelial Cells (HAEC)	C0065C
	Gibco Human Umbilical Vein Endothelial Cells (HUVEC)	C0035C
	Other primary cells and cell lines	Various
Instruments	Countess 3 Automated Cell Counter	AMQAX2000
	EVOS M3000 Imaging System	AMF3000
	Attune NxT Flow Cytometer, blue/red/violet/yellow	A29004

 Find out more about Nunc cultureware at
thermofisher.com/cultureware

thermo scientific