

Analyzers that meet your demands



Analyzers for Natural Gas and Natural Gas Liquid

Over 20 standard analyzers available that meet your demands - tested and ready to start-up at installation. These cover BTU determination, composition analysis, and impurities (e.g., sulfur species).

Single Channel Analyzers - Lowest Price

Single Channel systems analyze a single sample at a time. Analyzers can be gas only, liquid only, or gas or liquid injection.

Dual Channel Analyzers - Best Value

Dual Channel analyzers can be configured to provide high sample throughput (e.g., dual-sample analysis by GPA 2261) for greater productivity. Dual Channel systems are also available to provide more data per sample — the GPA 2286 system provides both BTU content and detailed hydrocarbon profile for a single sample/injection.

Super Combination Systems - Most Capability

For maximum flexibility in a single analyzer, Super Combination systems handle both NG and NGL samples to meet the requirements of all four GPA methods, including O_2/N_2 and He/H₂.

An additional channel for sulfur detection by ASTM D6228 can be added to any system.

	Thermo Fisher Scientific NG and NGL Analyzers		Natural Gas				Natural Gas Liquid	
Part number	Description	Samples	GPA 2261	O ₂ /N ₂	He/H ₂	GPA 2286	GPA 2177	GPA 2186
NATGAS010011	Natural Gas Analyzer with early back flush C6+ peak	1	~					
NATGAS010021	Natural Gas Analyzer with early back flush C6+ peak plus O2/N2	1	~	~				
NATGAS020271	Natural Gas Analyzer with early back flush C6+ peak and H2/He	1	~		~			
NATGAS020281	Natural Gas Analyzer with early back flush C6+ peak plus O2/N2 and H2/He	1	~	~	~			
NATGAS010031	Liquefied Natural Gas with early back flush C7+ peak	1					~	
NATGAS010211	Natural Gas with C6+ plus O2/N2 or LNG C7+ backflush	1	~				~	
NATGAS010041	Natural Gas extended hydrocarbon channel only	1				✓ *		
NATGAS010051	Liquified Nautural Gas extended hydrocarbon channel only	1						✓ *
NATGAS010221	Natural Gas or LNG Analyzer extended hydrocarbon channel only	1				✓ *		✓ *
NATGAS020231	Natural Gas with C6+ and extended hydrocarbon channel	1	~			~		
NATGAS020261	Natural Gas with C6+ plus O2/N2 and extended hydrocarbon	1	~	~		~		
NATGASC20211	Natural Gas or LNG with C6+ or C7+ and extended hydrocarbon	1	~			~	~	~
NATGASC20221	Natural Gas or LNG with C6+ plus H2/He or C7+ and extended hydrocarbon	1	~		~	~	~	~
NATGAS020011	Dual sample Natural Gas with early back flush C6+	2	~					
NATGAS020021	Dual Sample Nat Gas Analyzer with back flush peak plus O_2/N_2	2	~	~				
NATGAS020031	Dual sample LNG Analyzer with early C7+ back flush peak	2					~	
NATGAS020041	Dual sample Nat Gas Analyzer extended hydrocarbon channel only	2				✓ *		
NATGAS020051	Dual sample Liquified Natural Gas extended hydrocarbon channel only	2						✓ *
NATGAS020211	Dual Sample Nat Gas or LNG with early back flush peak plus 02/N2	2	~	~			~	
NATGAS020221	Dual sample Nat Gas or LNG Analyzer extended hydrocarbon channel only	2				✓ *		✓ *
NATGAS040011	Four sample Natural Gas with early back flush C6+ peak	4	~					

Custom Solutions

Do you have a specific requirement that is not covered by a standard analyzer? We will create a custom solution for you. Our team of experts will work with you to define your requirements to create a solution that meets your needs. Custom solution systems come with the same comprehensive documentation package as our standard analyzers.

Natural Gas Analysis Pressure-Innovation-Solution

Hundreds of thousands of dollars change hands based upon the results of your assay - *pressure*. A plant many miles away sets up production based on your analysis - *pressure*. New contracts, more samples, and new labs open, promising faster turnaround times - *pressure*. There is

Thermo Scientific[™] natural gas analyzers are built to meet your demands and to provide innovative solutions to relieve this constant pressure. Designed by natural gas experts with over 30 years of industry experience, our analyzers are configured and tested to the exact methods you follow:

constant pressure to be accurate, responsive, and competitive.

- GPA 2261, 2177, 2186, 2286
- ASTM D1945, D1946
- DIN 51872-5
- ISO 6974













It Starts with the Gas Chromatograph

The innovation in the TRACE[™] 1300 Series GC Gas Chromatograph will change the way you manage your GC systems - forever *and for the better*.

"Instant Connect" modularity, a patented Thermo Fisher Scientific technology, has completely revolutionized GC design. These modules incorporate all relevant pneumatic hardware and electronic parts necessary for making the injector or the detector a fully self-sufficient subunit of the instrument. The modules are plugged in to the top part of the GC, automatically configured into the system, and connected to the gas supply lines. Installing a module takes only two minutes.

- No special training required reduce training costs
- No need for a service call reduce cost of ownership
- No extended downtime higher uptime and more samples per month

Valve, Plumbing and Valve Ovens

Each natural gas analyzer is delivered ready to run your samples. Our solutions maximize performance and minimize footprint.

Compact Valve Oven - simplified configurations

- Valves: Two heated or unheated
- Maximum temperature: 160 °C
- Easy access bulkhead In/Out fittings


Auxiliary Oven - increased capacity, flexibility and capability -

- Valves: Up to six heated and two unheated rotary valves
- Maximim temperature: 350 °C
- Room for multiple columns
- Optional separate column oven module for additional isothermal zone
- Optional two detector module on top simultaneous and independent 4-channel capability


We deliver more to meet your demands


Thermo Fisher Scientific natural gas analyzers are designed, built and tested to meet method requirements, withstand the rigors of high throughput production labs, and increase analytical performance.

- Accurate and Reproducible Confidence in results every time
- Robust Built to optimize high sample loads
- Reliable Maximized throughput and minimized downtime

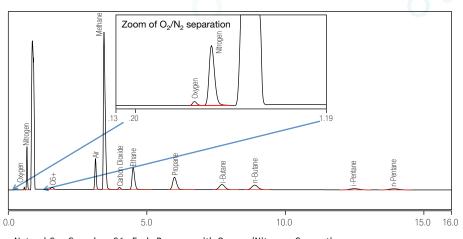



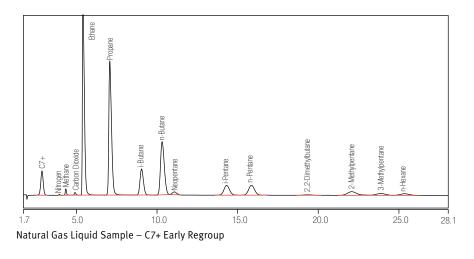

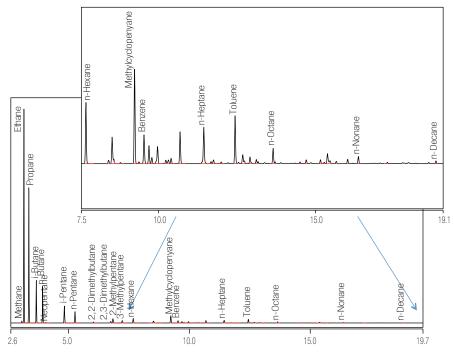

Extended Chromatogram for GPA 2186 (inset zoom of C6+)



©2016 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries. This information is presented as an example of the capabilities of Thermo Fisher Scientific products. It is not intended to encourage use of these products in any manners that might infringe the intellectual property rights of others. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

 Africa
 +43 1 333 50 34 0
 Denmark
 +45 70 23 62 60

 Australia
 +61 3 9757 4300
 Europe-Other
 +43 1 333 50 34 0

 Australia
 +61 3 9757 4300
 Finland
 +358 10 3292 200

 Belgium
 +32 53 73 42 41
 France
 +33 1 60 92 48 00

 Brazil
 +55 11 2730 3006
 Germany
 +49 6103 408 1014

 Canada
 +1 800 530 8447
 India
 +91 22 6742 9494

 China
 800 810 5118 (ree call domestic)
 400 650 5118
 +39 02 950 591

Japan +81 6 6885 1213 Korea +82 2 3420 8600 Latin America +1 561 688 8700 Middle East +43 1 333 50 34 0 Netherlands +31 76 579 55 55 New Zealand +64 9 980 6700 Norway +46 8 556 468 00

 Russia/CIS
 +43
 1
 333
 50
 34
 0

 Singapore
 +65
 6289
 1190
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90
 90</

