thermo scientific

Confident quantitation

A new level of confidence



Thermo Scientific TSQ Fortis triple quadrupole mass spectrometer




Fast and robust quantitation for every Environmental and Food Safety workflow




With evolving challenges in environmental and food safety applications, it's important to have the right platform for performing confident quantitation of multiple analytes in a complex mixture. The Thermo Scientific[™] TSQ Fortis[™] triple quadrupole mass spectrometer provides superior robustness—enabling round-the-clock operation and the ability to report reliable results, rapidly.

- Increased robustness—Active Ion Management Plus (AIM+) includes the Matrix Separator Ion Guide (MSIG)—ensures extended system uptime by keeping the ion path clean
- Maximum productivity—fast turnaround times while ensuring excellent quantitative results for extended lists of contaminants
- **Designed for more ease**—comprehensive workflows enabled by software that ensure confident analysis

Quantitation of contaminants in environmental samples at regulatory levels

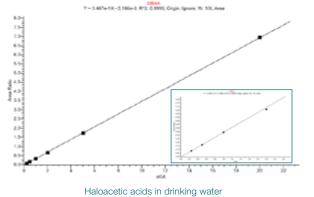

Quantitation of pesticides in food samples at regulatory levels

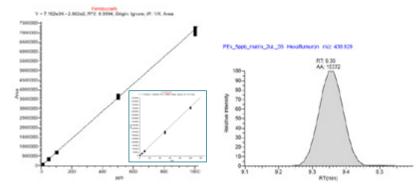

High throughput with ease of use for **Clinical Research** and **Forensic Toxicology**

Singlechannel





Clinical Research and Forensic Toxicology Workflow

Multichannel







The TSQ Fortis MS provides a seamless user experience for clinical research laboratories that focus on processing and analyzing a large number of samples. Designed for ease, the TSQ Fortis MS platform enables confident quantitation of hundreds of samples in complex matrices, every day.

- Enhanced robustness—AIM+ and the MSIG enable throughput by extending system uptime
- Usability-simplified user maintenance and intuitive software experience
- High throughput—analyze hundreds of samples and numerous compounds with short turnaround times


Consistent performance under challenging conditions

Precision, accuracy and sensitivity for demanding applications

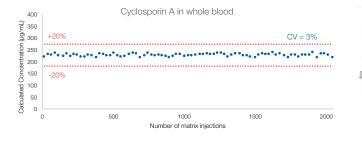

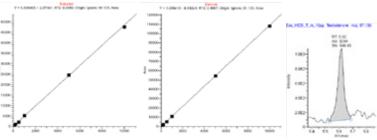


Reliability and sensitivity to meet **Pharmaceutical** requirements

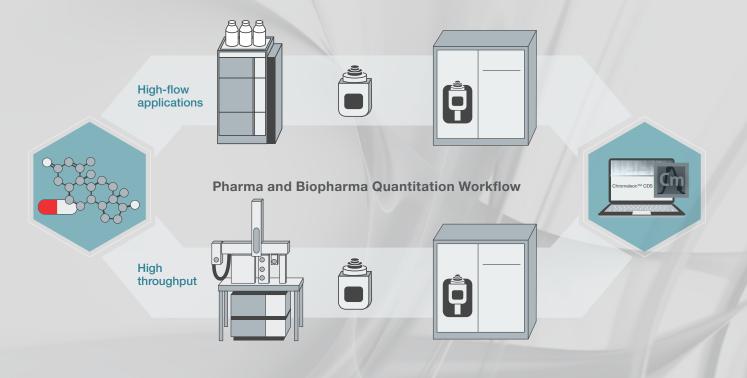

Pharmaceutical laboratories require ultimate reliability and the ability to achieve desired sensitivity to meet analytical demands. The TSQ Fortis MS provides more ease and comprehensive workflows empowering the daily quantitative requirements of pharmaceutical labs.

- **Outstanding reliability**—robust and reproducible data for all molecule types provided by AIM+ and the MSIG
- Sensitivity to meet demands—segmented quadrupoles facilitate ion transmission with system-to-system reproducibility
- High-throughput quantitation—multiple liquid chromatography (LC) options combined with the speed of the TSQ Fortis MS help analyze hundreds of samples every day

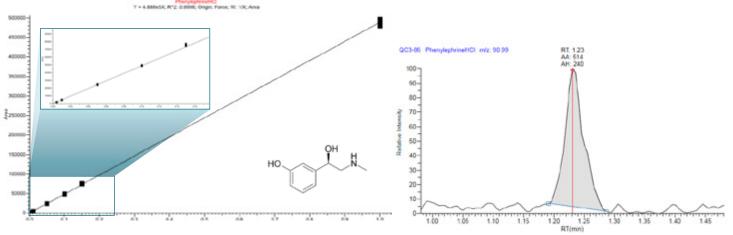

Ultimate productivity for QA/QC workflows

Companion solutions for the TSQ Fortis mass spectrometer

Workflow-driven software, solutions, and support streamline your applications

Integrated workflow-driven software and solutions streamline everyday tasks from instrument optimization and troubleshooting to routine method development and results reporting.



Thermo Scientific[™] TraceFinder[™] software

Comprehensive quantitative workflows for all applications from method development to report generation



AppsLab method library

A comprehensive repository for application methods, data sets, and application notes



Thermo Scientific[™] Chromeleon[™] software

Offers enterprise instrument control, data processing, and ability to address regulatory requirements

Our global Centers of Excellence, in-your-lab training using your methods, and AppsLab Library of Analytical Applications, ensure your success.


Chromatography solutions for MS

The TSQ Fortis mass spectrometer can be paired with a range of high-performance Thermo Scientific" inlet systems.

- Thermo Scientific[™] Vanquish[™] Flex, Horizon and Duo UHPLC systems
- Thermo Scientific[™] Transcend[™] II system with Multi-channel and TurboFlow[™] technology
- Thermo Scientific[™] Prelude SPLC[™] system with Multiplexing and TurboFlow[™] technology
- Thermo Scientific[™] Dionex[™] Integrion[™] HPIC[™] system



Experience confident quantitation: TSQ Fortis MS with **AIM+** Technology

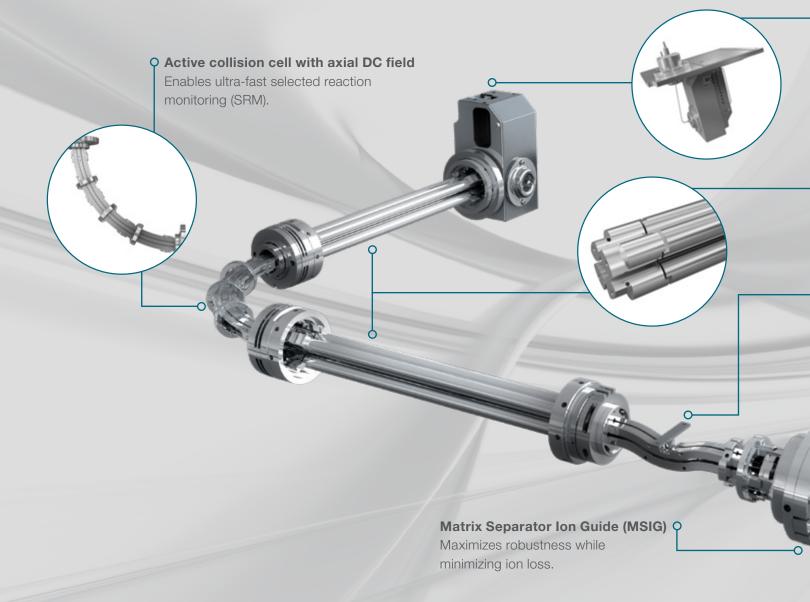

Robust, innovative design Benefits of the TSQ Fortis mass spectrometer

More samples, more compound classes

The reliable instrument design, intuitive software, and comprehensive workflow solutions developed with the TSQ Fortis triple quadrupole mass spectrometer help ensure superior productivity for quantifying hundreds of compounds of all types—in any matrix, by any user.

Robust quantitative productivity

The TSQ Fortis triple quadrupole mass spectrometer is designed to meet demands day after day, for every sample type in a variety of matrices. In addition to robustness and reproducibility, the TSQ Fortis mass spectrometer provides unprecedented ease of use. • Enhanced dual-mode discrete-dynode electron multiplier detector

Extends multiplier lifetime with increased surface area. Maintains excellent linearity and dynamic range.

• QR4 Segmented Quadrupole mass filter with hyperbolic surfaces

Simplifies day-to-day operation by ensuring consistent performance and reproducible results across instruments and over time for increased productivity.

O lon beam guide with neutral blocker

Provides robust performance and maximum sensitivity by blocking neutrals and efficiently transmitting ions.



• OptaMax NG ion source APCI ready

Robustness by design

Active Ion Management

The next step in precision design delivers the ultimate

Thermo Scientific[™] OptaMax[™] NG ion source housing

segmented quadrupoles and enhanced RF electronics

is designed with robustness in mind. Regardless of

your application, the TSQ Fortis mass spectrometer will deliver results reliably and consistently, day after day.

Plus (AIM+)

Automates all gas and voltage connections for ease of use. Provides optimizable spray position for ultimate performance in HESI or APCI mode.

P Ion transfer tube and sweep cone Provides maximum robustness and tool-free maintenance.

Confident data: technology translated to results

The TSQ Fortis mass spectrometer equipped with the neutral blocker, ion beam guide, and segmented quadrupoles is designed to deliver robust performance.

Speed with selectivity and sensitivity

The TSQ Fortis mass spectrometer delivers excellent quantitative performance in both selected reaction monitoring (SRM) and high-resolution SRM (H-SRM) modes, even when operated at short dwell times.

Your everyday partner for confident quantitation

The TSQ Fortis MS system software is an intuitive package designed to ensure ultimate productivity.

- Simple maintenance—easy clean-up of ion transfer tube without breaking vacuum
- High speed and automation—address hundreds of samples per day, or hundreds of transitions per sample
- Intuitive drag-and-drop method editor—software with application-based templates helps simplify method development and operation

Transforming triple quadrupole mass spectrometers

The TSQ Fortis mass spectrometer is one in a family of latest-generation Thermo Scientific[™] triple quadrupole mass spectrometers that offer confident quantitation every day, for every lab and every user.

These instruments are built on a foundation of shared, state-of-the-art hardware and software components. This makes it easier to transfer methods from one instrument to another when methods progress from single, sample-based experiments to validation and routine high-volume screening and quantitation.





Almanac web-based monitoring and management

The ability to remotely access and manage your LC-MS and GC-MS systems lets you focus on other important tasks during your day. Using the Almanac phone App, Instrument Connect or web-based application, you can check real-time system status and acquisition, set up automated e-mails to notify you of a completed acquisition or error, schedule instrument access, monitor utilization, or submit information to service personnel to aid in instrument diagnostics. thermofisher.com/almanac



Technical and online support: peak performance for your instruments

Helping you keep your instruments running at peak performance is our goal. Whether you're looking for an instrument manual or spare parts, want to submit a repair request, or check on the status of your warranty or service contract, we have every support option you're looking for. thermofisher.com/technicalresources

Protecting your investments: unparalleled laboratory services

Unity[™]Lab Services provides a single source for integrated lab service, support, and supply management. Our customized service offerings and world-class service experts have the flexibility and experience to effectively address your laboratory's business needs. We provide a complete portfolio of services and support solutions designed to help you improve productivity, reduce total cost of ownership, and ensure performance throughout your laboratory—from instrument and equipment acquisition to disposition. unitylabservices.com

Find out more at thermofisher.com/Fortis



For Research Use Only. Not for use in diagnostic procedures. © 2018 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. **BR65172-EN 0819M**