



Centrifugation

# Secured at every turn


## New SECURE-Spin mode for the General Purpose Pro Centrifuge Series\*

### Add advanced and secure control features to your centrifuges

In regulated environments, it is crucial to assign user rights and maintain data management in a 21 CFR Part 11-compliant manner. This ensures traceable documentation to help maintain sample and process safety, to provide audit trails, and to record data securely. The new Thermo Scientific™ SECURE-Spin mode takes the Thermo Scientific™ General Purpose Pro Centrifuge Series\* to the next level with accurate data collection, reporting, and archiving, as well as management and protection of process data. SECURE-Spin mode adds advanced functionality to your centrifuge, allowing users to:

- Create up to 120 users registered with unique usernames and passwords
- Assign individual privileges and roles to each user with different authorities
- Track all user changes to programs or settings with the audit trail

### SECURE-Spin mode supports 21 CFR Part 11 compliance

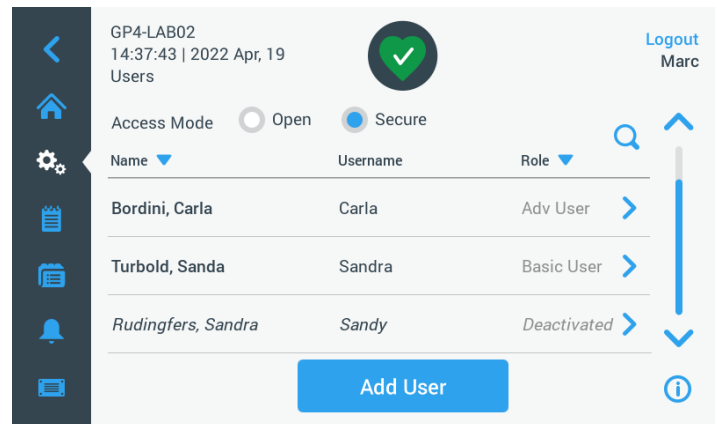

\* Applies to Thermo Scientific™ Multifuge™ X Pro and Thermo Scientific™ Sorvall™ X Pro centrifuges only.

## Extending Centri-Log Plus software connectivity to General Purpose Pro centrifuges

Functions of centralized monitoring of an unlimited number of centrifuges with accurate audit trails using Thermo Scientific™ Centri-Log™ Plus Data Management Software

- Identify each centrifuge by a serial number, assigned unit name, or IP address on the same network
- Gather real-time run data, including speed, time, temperature, alarm messages, and more, displayed on one PC, for all networked centrifuges
- Automatically record run data that will be stored directly on the assigned PC or network drive in a .csv file format
- Retain any changes in centrifuge, user, and program settings with additional timestamped information on the user who created, modified, or deleted the record



### Ordering information

Description	Cat. No.
<b>Software and accessories</b>	
SECURE-Spin mode	75009990
Centri-Log Plus Software** (requires SECURE-Spin mode)	75007742
Sample tracking kit* (requires Centri-Log Plus software)	<ul style="list-style-type: none"> <li>• For floor-standing centrifuges (barcode scanner, holder, and cable) 75007740</li> <li>• For benchtop centrifuges (barcode scanner and cable) 75007743</li> </ul>
Network access kit* (cable)	75007741

\* Product sold separately.

† Centri-Log Plus software version must be at least 2.12.0.



Learn more at [thermofisher.com/securespin](https://thermofisher.com/securespin)