

August 2018

Subject: Change notification for ClipTip 1250 µl pipette tips

Dear Customer,

Thermo Fisher Scientific is pleased to announce that we have made enhancements based on customer feedback to the Thermo Scientific™ ClipTip™ 1250 µl Pipette Tips to improve performance of the tips when used with E1-ClipTip pipettes.


The improvements include:

- Thinner and pointier shape that allows more precise aspiration from the bottom of vessels
- 5 mm longer tip, enabling better reach
- New filter raw material and a shorter filter that improves stepping performance

The changes are implemented on the below product codes:

Package description	Cat. no.
ClipTip 1250, Racked, 96 tips/rack, 8 racks/pack	94410810
ClipTip 1250, Racked, Sterile, 96 tips/rack, 8 racks/pack	94410813
ClipTip 1250, Reload Insert, 96 tips/insert, 8 inserts/pack	94410817
ClipTip 1250 Reload Insert, Sterile, 96 tips/insert, 8 inserts/pack	94410818
ClipTip 1250, Racked, Filtered, Sterile, 96 tips/rack, 8 racks/pack	94420813
ClipTip 1250, Reload Insert, Filtered, Sterile, 96 tips/insert, 8 inserts/pack	94420818

OLD CT1250 NEW CT1250

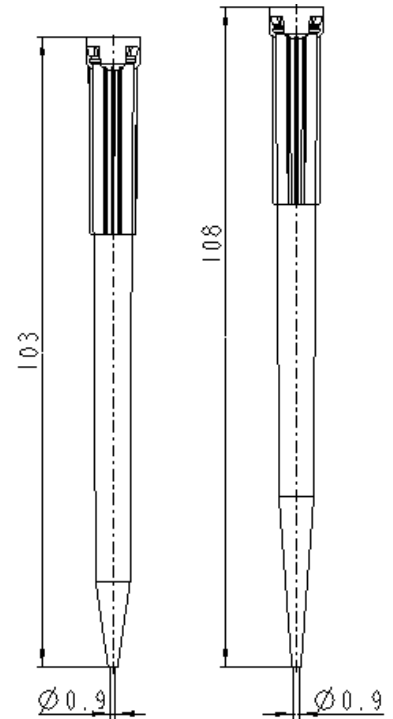

Changes were implemented from lot number 18205xx forward. While the product codes stay the same, we have included a sticker on the product packages to indicate the design change. Additionally, there is a package insert card inside the sales units of the new format to explain the design change.

There are no changes to product specifications (accuracy and precision values with E1-ClipTip pipettes), descriptions or pricing. Filter raw material changed only on physical properties to a material that has been used for years in our other pipette tips.



It's recommended that the E1-ClipTip pipette calibration is checked when switching from old tips to new tips. However, based on our internal testing, it's likely that no adjustments should be needed.

The old format of the tip will be no longer produced. It will be available as long as stock remains and then we will move forward to the new design.

If you have any questions or require additional information, please contact our product support: info.pipettes@thermofisher.com.

Thank you for your continued support of Thermo Scientific products and services.

Sincerely,

Maiju Palokangas
Senior Product Manager
ClipTip products and Core Pipettes
Liquid Transfer

Marina Elf
Quality and Compliance Manager
Laboratory Products Division

