

# ToGo Viewer Software

Publish and share your 3D data analysis results

You and your organization are performing acquisitions and analysis on a regular basis with our Amira Software and Avizo Software products.


You need to be able to share the results of data analysis with your customers, your business partners, or within your organization.

Some aspects of such analysis may not always be easy to understand using only pictures, charts, spreadsheets or even videos.

In order to fully understand the results of the analysis, nothing is better than an interactive virtual simplified representation of the processed data.

## What is ToGo Viewer Software

Using the ToGo Publisher, which is part of Amira and Avizo 3D Pro, Avizo 3D for Industrial Inspection and Amira 3D for Cell Biology, you can publish the results of your analysis using a special Avizo data format, securely embedding all your project data and information in a single file. This ToGo file format includes the dataset itself, its properties, animations, colormap editor and tables.

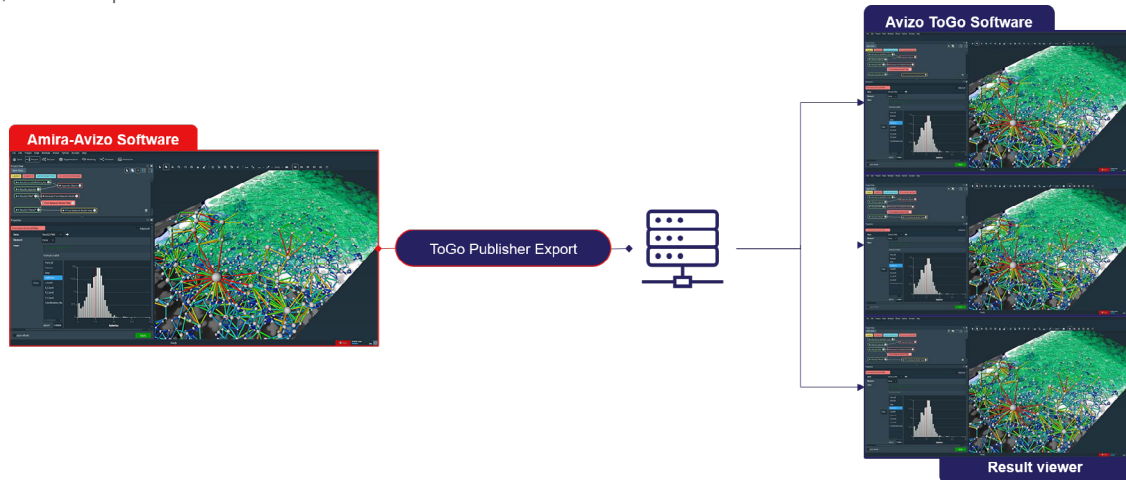

This special data file can then be loaded by ToGo viewer applications, on any computer by your customers / partners / collaborators.

ToGo Viewer software is a limited version of Amira and Avizo software, which can only load ToGo project files, and which provides interactive visualization features, with direct manipulation in the viewer window or by controlling the colormap or editing the exposed properties of display modules. ToGo Viewer software can also play animated movies generated with the project and embedded into the ToGo file.

## Providing ToGo Viewer Software to your users

Simply send the download link for the ToGo Viewer Software installer to your data users, and they will be able to access the published analysis results with ease:

<https://www.thermofisher.com/software-em-3d-vis/customerportal/download-center/togo-viewer-installers/>

## Help and Support

- Please make sure that your ToGo Viewer users have a sufficient hardware configuration to run the product and potentially large datasets, i.e. a modern graphic board and enough memory.
- Please refer to the ToGo Publisher chapter from the Amira or Avizo Software Help Menu for more details.
- ToGo Viewer Software users have access to the Amira and Avizo Software documentation, with specific introduction to ToGo Viewer Software.
- Should you or your user have questions regarding the use of ToGo Viewer software, please contact us through: <https://msdcss.powerappsportals.com/>

Find out more at [thermofisher.com/amira-avizo](https://www.thermofisher.com/amira-avizo)

**ThermoFisher**  
S C I E N T I F I C